


INSTRUCTIONS FOR USE


TWO-HANDED APPLICATION

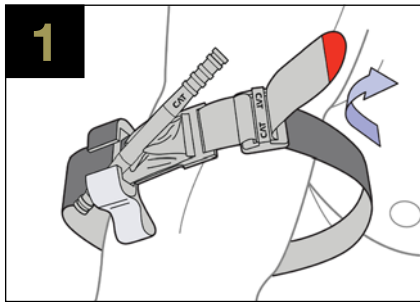


Route the band around the limb, pass the red tip through the slit of the buckle, and position the **CAT** 2-3 inches above the bleeding site directly to the skin.*

OR


ONE-HANDED APPLICATION

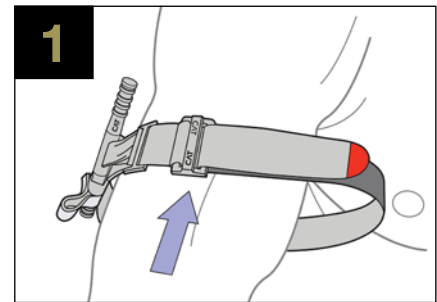


Insert the injured limb through the loop in the band and position the **CAT** 2-3 inches above the bleeding site directly to the skin.*

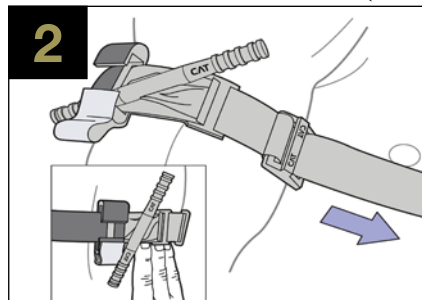


Pull band tightly and fasten it back on itself all the way around the limb, but not over the rod clips. Band should be tight enough that tips of three (3) fingers cannot be slid between the band and the limb. If the tips of three (3) fingers slide under band, retighten and re-secure.

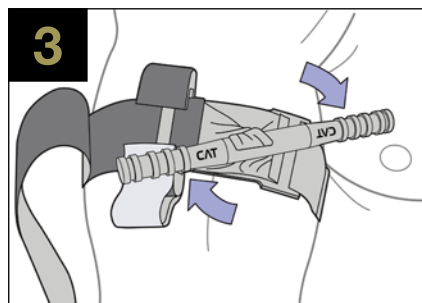

Twist the rod until bleeding has stopped.

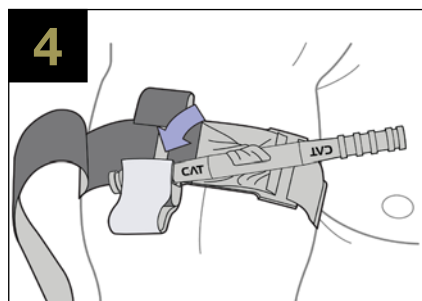

Secure the rod inside a clip to lock it in place. Check for bleeding and distal pulse. If bleeding is not controlled, or distal pulse is present, consider additional tightening or applying a second **CAT** above and side-by-side to the first. Reassess.

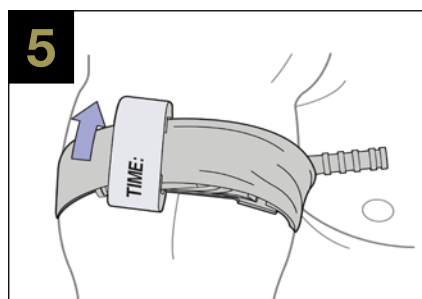


Route the band between the clips and over the rod. Secure rod and band with TIME strap. Record time of application.


*If you cannot be sure or cannot take the additional time to examine where the bleeding is coming from based on the situation, the **CAT** can be effectively applied over clothing as high on the arm or leg as possible. The **CAT** must NOT be applied over solid objects within the clothing. As soon as the situation permits, the injured limb should be evaluated and the **CAT** re-positioned 2-3 inches above the injury directly to the skin.




STORING THE CAT IN THE QUICK LAUNCH CONFIGURATION

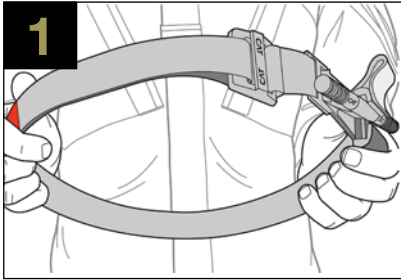

Pass the tip through the slit in the buckle. Pull 8 inches of band through, fold it back and adhere the band to itself.

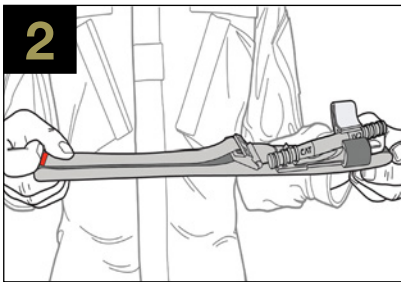

Flatten the loop formed by the band. Place the buckle in the middle of the flattened band.

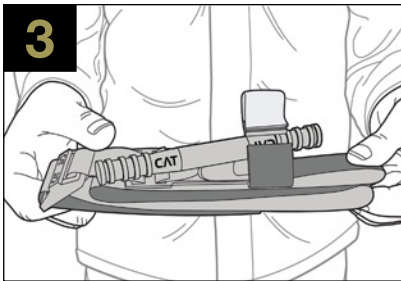

The CAT is now ready to be placed in your medical kit.

STANDARD USE GUIDELINES

Intended Use of the Combat Application Tourniquet® (CAT) is to occlude blood flow of an extremity in the event of a traumatic wound with significant hemorrhage. Clinical benefit of the CAT is the control of a life-threatening bleed which results in less blood loss. Intended user of the CAT should only be used as directed by user's military service component guidelines, EMS authority, or under the supervision of a physician, and in strict compliance with these Instructions for Use. Target population for the CAT is anyone with major bleeding to the upper or lower extremities, which is life threatening and cannot be stopped otherwise. The CAT has a Maximum Operating Circumference of 35in (88.9cm). The use of any tourniquet for longer than 2 hours may lead to permanent neurological or muscular damage. If you cannot be sure or cannot take the additional time to examine where the bleeding is coming from based on the situation, the CAT can be effectively applied over clothing as high on the arm or leg as possible. The CAT must NOT be applied over solid objects within the clothing. As soon as the situation permits, the injured limb should be evaluated and the CAT re-positioned 2-3 inches above the injury directly to the skin. The CAT should be stored in a cool, dry place away from direct sunlight. Re-use of this single use device may introduce risk of infection. Dispose of used CAT in accordance with local regulations for biomedical waste. Please report all serious incidents in relation to the CAT to the manufacturer and appropriate regulatory authority.

Manufacturer	Importer	Medical device	Medical device	Consult Instructions for Use	Date of Manufacture
Authorized representative in the European Union	Authorized representative in Switzerland	Authorized representative in the United Kingdom			
Keep dry	Keep away from sunlight	Non-sterile	Unique Device Identifier	Lot Number	

Combat Application Tourniquet® and CAT are trademarks registered in the United States and certain other countries.

U.S. Patent Nos. 7,842,067; 7,892,253;
8,888,807; 10,016,203;
Canada Patent No. 2,569,550;
EU Patent Nos. 1,753,344; 2,425,786;
Israel Patent No. 179769;
Hong Kong Patent No. 1,163,479;
Other Patents Pending

Manufactured by C-AT Resources, LLC
4203 International Dr., Rock Hill, SC 29732
United States of America
803.325.9300




NSN 6515-01-521-7976



REF 30-0001

Tactical Black

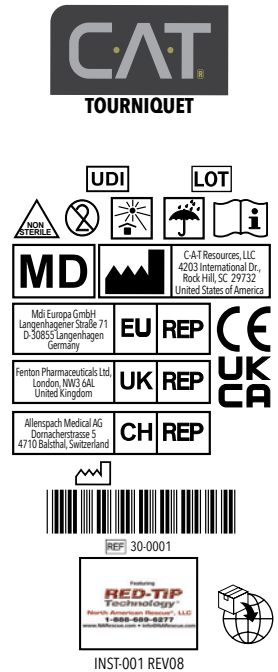

Exclusive distribution in the U.S. by:
North American Rescue®, LLC
1-888-689-6277
www.NARescue.com
info@NARescue.com



TACTICAL BLACK
NSN 6515-01-521-7976



TRAINER BLUE
NSN 6910-01-560-2972



EMS ORANGE
NSN 6515-01-612-7129