

# ***Ballistic Protection***



BALLISTIC PROTECTION | BALLISTIC VEST SYSTEMS | MEDICAL SUPPLIES



# OVERVIEW | Body Armor


The term body armor is generally used to describe vests designed to provide ballistic protection to the vital organs in the torso, armor has typically been associated with tactical medics and rescue task force members whose primary objective is to provide medical care in high threat environments. With the increase of active shooter and mass casualty incidents in today's field of operation, Fire & EMS personnel are realizing the critical need for ballistic protection on a daily basis.

*Though Fire & EMS personnel have been limited to adapting body armor traditionally designed for the military & law enforcement, North American Rescue is now confronting this challenge by designing ballistic solutions specifically for the unique needs of the Fire & EMS community.*

## Soft & Hard | 2 Main Types of Body Armor

**SOFT ARMOR** offers critical protection against handgun threats along with the capability for extended daily wear. NAR provides an optimal over-the-uniform soft armor vest that features flexible layers of ballistic resistant materials constructed from a synthetic Aramid fiber with a very high strength to weight ratio. These super strong fibers are woven into multiple layers that are stacked to create an armor panel. The number and orientation of each layer influences its overall performance. When a handgun bullet strikes the armor panel, it gets caught in a web of high performance fibers that stretch and disperse the impact energy causing the bullet to deform or "mushroom". This dispersion of energy not only prevents bullet perforation, but also helps minimize the possible injury to underlying tissue.

How Soft  
Armor  
disperses  
energy  
from a  
handgun  
projectile...

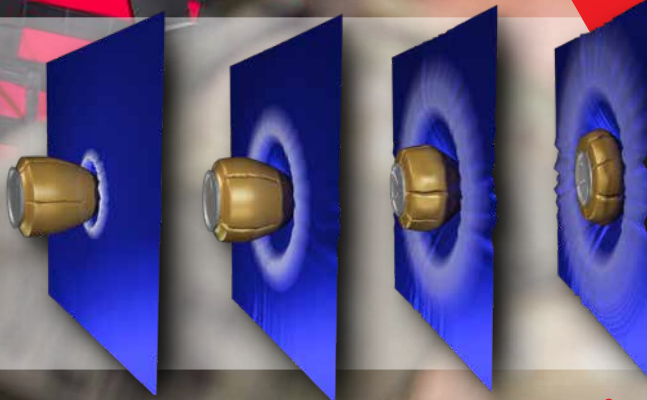

**HARD ARMOR** is designed to provide a higher degree of protection from rifle threats. Constructed from a variety of hard ballistic-resistant materials such as ceramic, compressed UHMWPE (ultra-high-molecular-weight-polyethylene), metallic plates, or a composite of multiple materials, each of these rigid materials offer their own unique properties. Ceramic materials, for example, are required for so called "armor piercing" bullets or special threats. Hard armor plates are classified as either Stand-Alone (no soft armor required) or ICW - in conjunction with (soft armor backing panel required) and rated as single hit or multi-hit. ▶

## NO armor is 100% bullet proof!

One truth stands out when it comes to protective gear – ballistic protection is NOT created equal. Differing armors are designed to protect against specific types of ammunition, and no body armor will stop every threat. So when choosing a ballistic solution, consider the threats most likely to be faced in your area of operation. At a bare minimum, armor should protect against handguns.



# NIJ 0101.06 Ballistic Standards

## Special Threats

In addition, the FBI and DEA have developed their own stringent protocols to test body armor that exceeds the NIJ standards.

### 5.56mm ammunition (M855)

was designed for military use and has a steel core. It is the most common and widely used AR ammo on the market today. To provide requisite protection, a ceramic component must be included in the hard armor plate.



## Choose the right ARMOR for YOU...

There are several key considerations in choosing the right body armor: level of protection, weight, comfort, mobility, and cost. It is important to choose an armor system that will provide protection but also enhance your ability to perform your job. Particularly for Fire & EMS personnel, it is critical that the system facilitates your ability to provide the life saving skills and equipment.

*North American Rescue has the optimal ballistic protective solution designed as a system specifically for Fire & EMS personnel and the duties they must perform.*

The **National Institute of Justice (NIJ)** is the government agency responsible for establishing performance standards and testing certification for body armor. These body armor classifications: levels IIA, II, IIIA, III and IV correspond with the weaponry they protect against. NIJ certified armor ensures ballistic performance and the completion of proper testing standards.

The NIJ conducts two types of ballistic testing: Perforation/BackFace Signature (P-BFS) and Ballistic Limit (V50). No complete or partial perforation of the armor panel is allowed during testing. In addition, no more than 44mm (1.73 inches) of indentation is allowed into the ballistic clay backing utilized during shoot testing. The back face signature translates to the upper tolerable limits of blunt force injury to the user. The ballistic limit (V50) is the velocity at which the projectile is expected to perforate the armor 50% of the time. It is important to ensure that any armor package has been certified by an NIJ approved laboratory.

## Level IIA:

Tested to stop 9mm and .40 S&W ammunition fired from short barrel handguns

No rifle ammunition protection

## Level II:

Tested to stop 9mm and .357 Magnum ammunition fired from short barrel handguns

No rifle ammunition protection.

## Level IIIA:

Tested to stop .357 SIG and .44 Magnum ammunition fired from longer barrel handguns.

No rifle ammunition protection

## Level III:

Tested to stop 7.62mm FMJ (full metal jacket) lead core rifle ammunition

## Level IV:

Tested to stop .30 cal steel core armor piercing rifle ammunition

### NIJ Threat Projectiles & Related Armor Classifications





## MKIII TACTICAL RESPONDER BALLISTIC VEST *with Cummerbund*

- Developed for providing Fire and EMS personnel in active shooter and mass casualty incidents with Level IIIA Soft Armor Ballistic Protection
- Allows for full range of motion while providing superb coverage area
- Includes front and rear NIJ 0101.06 Level IIIA ballistic panels and two integrated plate pockets designed to accommodate additional 10 in. x 12 in. or 11 in. x 14 in. hard armor plates (each sold separately)
- Featuring extra side protection with Level IIIA armor cummerbund equipped with easy hook and loop attachment to the vest for rapid user adjustment
- MOLLE-style PALS webbing on the front, back and around the waist allows for attaching extra gear and for unit customization
- Rapidly adjustable armor system customizes to any team member's frame within seconds
- Can be stored in and deployed from all standard emergency vehicles

### MKIII TACTICAL RESPONDER VEST with CUMMERBUND

COLOR	ITEM #
Black	70-0251
Blue	70-0253
Red	70-0254
Coyote	70-0252

#### DIMENSIONS

Vest: H 19 in. x W 12.5 in. x D 4 in.

Vest with Cummerbund:

H 19 in. x W 15.5 in. x D 4 in.

Maximum Waist Size: 62 in.

Total Weight: 9 lb 6 oz

#### BALLISTIC SPECIFICATIONS

Protection Type/Threat Level: NIJ 0101.06 IIIA  
Special Threat Tested

Construction: 3000D Twaron® Microlam |  
Twaron® 2-Ply Unidirectional Fabric

W.R.A.P.™ - Water Repellant Armor Perimeter  
(patent pending)

Protective Cover: Brookwood 70D  
with Laminated Polyurethane

Thickness: 0.275 in. | AD= 1.17 psf

Model#: Fusion IIIA AFC-T3A-525  
(NIJ Compliance Product List)

Testing Certification: NIJ 0101.06  
and Special Threat Tested

#### Special Threats:

.44 Mag V50: 1615 fps

.357 Sig V50: 1735 fps

FN SS197 40 grain VO

Win SXT 9mm+P+ VO

7.62 x 25 Tokarev 85 grain VO

.357 Sig WIN SXT 125 grain VO

Speer GD .357 Mag VO

2 grain RCC V50: 3013 fps

4 grain RCC V50: 2545 fps

16 grain RCC V50: 2225 fps

17 grain RCC V50: 2077 fps



MKIII

## RTF TACTICAL RESPONDER BALLISTIC VEST SYSTEM

*The practical solution for Fire and EMS personnel in active shooter and mass casualty incidents, with an emphasis on integrating Tactical Emergency Casualty Care...*

- NIJ 0101.06 Level IIIA Soft Armor Ballistic Protection equipped with medical supplies that address the leading causes of preventable death in the tactical environment
- Front and back hard armor plate pockets able to accommodate additional 10 in. x 12 in. or 11 in. x 14 in. hard armor plates (each sold separately)
- Includes Cummerbund-style Side Armor Protection with easy hook and loop attachment to the vest for rapid user adjustment - expandable up to 62 in.
- MKIII RTF Vest Tactical Responder Vest System includes:
  - 1 x MKIII Tactical Responder Ballistic Vest with Cummerbund
  - 1 x RTF Chest Pouch Medical Kit with Luminous Med ID Patch
  - 2 x Nylon Tourniquet Holder with C-A-T®
- ID patches are also available for separate purchase and are customizable upon special request.



### MKIII RTF TACTICAL RESPONDER BALLISTIC VEST SYSTEM

VEST	ACCESSORIES	ITEM #
Black	Black	80-1092
Blue	Black	80-1093
Red	Red	80-1095
Coyote	Coyote	80-1094

#### DIMENSIONS

Vest, Cummerbund & Medical Kit:  
H 19 in. x W 15.5 in. x D 7 in.

Maximum Waist 62 in.

Total Weight: 12.4 lb

Medical Kit (separately):  
H 7.75 in. x W 9.75 in. x D 4.5 in.  
Weight: 2 lb 10.3 oz

#### RTF CHEST POUCH MEDICAL KIT (also available separately)

COLOR	ITEM #
Black	80-0407
Coyote	80-1109
Red	80-0406

#### on MKIII CUMMERBUND:



#### RTF CHEST POUCH MEDICAL KIT Contains:





## Ballistic Rifle Plates

*Comfortable & lightweight...*

*The clear choice for an economical, multi-strike*

*NIJ certified advanced threat armor plate under 5 lb...*

### NIJ Level III Plus Enhanced **M7X**

- Protection Type/Threat Level: NIJ III Plus
- Configuration: Stand-Alone
- Construction: Ceramic /Polyethylene /Composite
- Protective Cover: Polyurethane
- Hit Rating: Multi-Strike
- Model#: DM855 (NIJ Compliance Product List)
- Positive Buoyant
- 7 Year Warranty



- Testing Certification: NIJ 0101.06 & Special Threat Tested | Tested to DEA standards against M855\* and 7.62 x 39 MSC

#### **Special Threats:**

- 5.56mm x 45mm 55gr M193 FMJ (NATO)
- 5.56mm x 45mm 62gr Ball (M855)\*
- 7.62mm x 39mm 123gr PS Ball (Mild Steel Core)
- 7.62mm x 51mm 147gr M80 FMJ



***\*Tested to stop M855 5.56 "Green Tip" - the cheapest, most widely available AR15 ammunition***



#### **M7X NIJ LEVEL III+ ENHANCED**

STYLE	WEIGHT	ITEM #
Full Size	5.15 lb	70-0255
Shooter	4.65 lb	70-0256

#### **SPECIFICATIONS**

W 10 in x H 12 in x D 125 in

## Ballistic Helmets



### Batlskin® Full Cut **VIPER A3 HELMET** with Modular Suspension System

- The cost effective ballistic and head impact protection solution
- Exceeds U.S. Army ACH impact standards and tested with multiple hits and across all temperature ranges
- Stops NIJ IIIA level threats (9mm and 44 mag) (NIJ 0101.06)
- Offers unsurpassed impact protection and ballistic rated hardware tested for ballistic performance from -40° F to 158° F
- Compatible with NVGs, Chem/Bio Masks & communication devices
- **Modular Suspension System (MSS)** features 4-point integrated harness with unique buckle system (no loose straps and Dial Ratchet Mechanism to adjust the fit band with one hand
- Adjustable nape pad ensures helmet stability
- Adjustable crown net allows sweat wicking, comfort and low pressure wear, providing a free flow of air to aid in cooling
- Low profile impact material does not contact the head in normal operation


#### **SPECIFICATIONS**

Weight: 3 lb

#### **VIPER A3 with MSS**

COLOR	SIZE	ITEM #	COLOR	SIZE	ITEM #
Black	S	70-0081	Coyote	S	70-0247
Black	M	70-0082	Coyote	M	70-0248
Black	L	70-0083	Coyote	L	70-0249
Black	XL	70-0084	Coyote	XL	70-0250



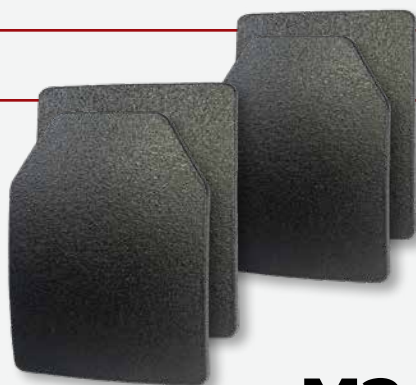

## M3 NIJ Level III Plus

*The lightest armor in its category!*

### M3 NIJ LEVEL III+

STYLE	SIZE	WGT	ITEM #
Full Size	10 x 12	3.15 lb	70-0257
Shooter	10 x 12	2.9 lb	70-0259

STYLE	SIZE	WGT	ITEM #
Full Size	11 x 14	3.97 lb	70-0258
Shooter	11 x 14	3.5 lb	70-0260

- Protection Type/  
Threat Level: NIJ III Plus
- Configuration: Stand-Alone
- Construction: Polyethylene
- Protective Cover:  
Polyurethane
- Hit Rating: Multi-Strike
- Model#: DKX MAX III  
(NIJ Compliant Product List)
- Positive Buoyant
- 7 Year Warranty
- Testing Certification:  
NIJ 0101.06  
& Special Threat Tested

#### Special Threats:

- .223 Rem 60gr Partition
- 7.62mm x 39mm 123gr  
Mild Steel Core (MSC)
- 7.62mm x 51mm 147gr  
M80 FMJ, V50 = 3288 fps
- 5.56mm x 45mm 55gr  
M193 FMJ (NATO)

## Ballistic Eye Protection

### NAR IPRO® Tactical Goggle System

- Fully enclosed, high-impact eye protection  
compatible with helmet and night-vision goggles
- 2 clear & 2 tinted lenses (neutral gray) with  
100% UV shield
- Hexo-vents for fog-free viewing
- Certified to ANSI Z87.1-2003 & EN 166:2001
- Includes cleaning kit & padded zipper case

#### SPECIFICATIONS

Case: L 2.6 in. x W 7 in. x D 3 in.  
Weight: 8 oz

#### NAR IPRO

##### ITEM #

ZZ-0109



## Nylon Accessories

### Rescue Task Force ID PATCH

Nylon Patches with Hook Backing  
for agency ID on Ballistic Vests

**RESCUE  
TASK FORCE**

**RESCUE TASK FORCE**

*Custom patches also available.  
Call Customer Service for more details.*

#### RESCUE TASK FORCE PATCHES

COLOR	SIZE	ITEM #
Red	3 in. x 11 in.	70-0200
Coyote	3 in. x 11 in.	70-0199
Black	3 in. x 11 in.	70-0198
Red	2 in. x 11 in.	ZZ-0537
Coyote	2 in. x 11 in.	ZZ-0510
Black	2 in. x 11 in.	ZZ-0509

### Nylon C-A-T® HOLDER

#### SPECIFICATIONS

L 6.25 in. x W 1.75 in. x D 1.5 in.  
Weight: 1.26 oz

#### C-A-T HOLDER

COLOR	ITEM #
Black	30-0057
Coyote	30-0044
Red	30-0128



### Nylon GEAR BAG

Single compartment bag for  
maximum capacity - large enough to  
store essential ballistic gear including  
vest & helmet



#### SPECS

L 30 in. x W 14 in.  
x D 14 in.  
Weight: 1 lb 4 oz

#### NYLON GEAR BAG

COLOR	ITEM #
Black	80-0547

### NAR SOFT LITTER HOLDER

#### SPECIFICATIONS

H 14 in. x W 14 in. x D 0.25 in.  
Weight: 9 oz

#### LITTER HOLDER

COLOR	ITEM #
Black	80-0556
Blue	80-0557
Red	80-0558



- Carries (2) NAR QuikLitters or (1) NAR QuikLitter XL (each sold separately)
- Allows deployment of litter from the bottom or top of the carrier
- Featuring a MOLLE-style attachment that connects to the rear of protective vest leaving front of vest open for other mission critical equipment





## INDIVIDUAL THROW KITS

- Compact vacuum-sealed pouch with essential point of wounding equipment provided for rapid deployment to on-scene responders during a Mass Casualty Incident
- Designed to treat bleeding from penetrating and other traumatic injuries
- Small cube space fits into a BDU pocket, vest pouch, or individual bag


### SPECIFICATIONS

H 7 in. x W 5 in. x D 2.5 in.  
Weight: 8 oz

### INDIVIDUAL THROW KITS

OPTIONS	ITEM #
Basic	80-1104
with Combat Gauze LE	80-1105

### INDIVIDUAL THROW KITS:

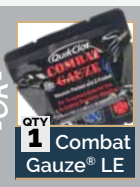

## MCI RAPID RESPONSE THROW KITS

*Hands Free, Point of Wounding Solution with the enhanced ability of treating a large number of injured patients in a large mass casualty incident*

- Rapid response sling bag designed to take the guesswork out of the management of a critical incident by providing the rapid distribution of necessary life saving equipment over a large scene
- Includes 6 Individual Throw Kits with essential life saving bleeding control equipment designed for rapid deployment to on-scene responders
- Nylon sling bag features exterior MOLLE-style attachments so the user can customize the bag according to their needs
- Provides ample storage for extra gear such as a Rescue Harness for casualty evacuation (sold separately)

### SPECIFICATIONS

H 12 in. x W 12 in. x D 6 in.  
Shoulder strap: 28 in. to 50 in.  
Weight: 6 lb 4.5 oz

### MCI RAPID RESPONSE THROW KITS:

OPTIONS	ITEM #
Basic	80-1106
with Combat Gauze LE	80-1107



## HASTY HARNESS

*Multi-purpose Harness*

- Indispensable gear for every webbing equipment toolbox
- Resourceful, multi-purpose harness for unlimited applications & maximum versatility
- Constructed of 22 ft loop of 1 in. Tech Tape (MBS: 4,400 lb) with opposing color thread to aid in visual safety inspection

### SPECIFICATIONS

Loop: L 10 ft 11 in. x W 1 in.  
Circumference - 21 ft 10 in.  
Min. Break Str: Tech Tape (4,400 lb)  
Sewn Sling (6,000 lb)

### HASTY HARNESS

COLOR	ITEM #
Red	90-0044
Black	90-0001

## EX-TRACT STRAP

*Point of Wounding Extraction*

- Point of wounding extraction
- Back pack carry
- Hasty harness for traditional and confined space lowering and hauling



### SPECS

Packaged:  
H 18 in.  
x W 3.5 in.  
x D 1 in.

Extended:  
L 58 in.

Weight: 9.3 oz

### EX-TRACT STRAP

ITEM #
90-0082



## Mass Casualty Incident **MCI WALK**

- All-in-one trauma treatment & evacuation platform for treatment of multiple trauma patients
- Choose from two different litter configurations: 8 non-rigid, disposable QuikLitters or 1 full-size, quad-fold, lightweight MedEVAC Litter
- Pocket system holds self- or buddy-aid point-of-wounding medical supplies
- Bag features back pack straps for easy transport from the vehicle to the point of need
- Ideal for use as a vehicle trauma kit

(80-0565/80-0458) Red MCI WALKs contain orange C-A-T® Tourniquets

(#80-0458) Red MCI WALK with MedEVAC Litter contains an orange MedEVAC Litter

### MCI WALK

COLOR	LITTER	ITEM #
Black	8 QuikLitters	80-0564
Red	8 QuikLitters	80-0565
Black	MedEVAC Litter (Black)	80-0457
Red	MedEVAC Litter (Orange)	80-0458

### SPECIFICATIONS

H 22.75 in. x W 16.5 in. x D 11.5 in.

Weight with 8 QuikLitters: 18 lb 12 oz

Weight with MedEVAC Litter: 29 lb 15 oz

### MCI WALK with QUIKLITTERS:

 QTY 8 C-A-T® Tourniquet	 QTY 8 ETD™ 6 in. Dressing	 QTY 5 Combat Gauze® LE	 QTY 8 Gauze, Compressed	 QTY 5 28F Naso Airway & Lube
 QTY 8 HyFin® Vent. Twin Pk	 QTY 8 NAR QuikLitter	 PAIR 8 Nitrile Gloves, Lg	 QTY 20 Combat Casualty Card	 QTY 5 Permanent Marker, Large

### MCI WALK with MEDEVAC LITTER:

 QTY 8 C-A-T® Tourniquet	 QTY 8 ETD™ 6 in. Dressing	 QTY 5 Combat Gauze® LE	 QTY 8 Gauze, Compressed	 QTY 5 28F Naso Airway & Lube
 QTY 8 HyFin® Vent. Twin Pk	 QTY 1 Rigid Litter MedEVAC	 PAIR 8 Nitrile Gloves, Lg	 QTY 20 Combat Casualty Card	 QTY 5 Permanent Marker, Large



Bag only option also available at [NARescue.com](http://NARescue.com)



online video available

< scan to view

## C-A-T® Combat Application Tourniquet Gen 7

- **Single route buckle** decreases application time to help decrease blood loss, allows for slack removal for proper TQ application
- **Applied same way** to both upper and lower extremities
- True **one-handed tourniquet**
- **Scientifically proven** 100% effective in occluding blood flow in both upper and lower extremities by the USAISR and recommended by CoTCCC



### SPECIFICATIONS

Packaged: H 6.5 in. x W 2.4 in. x D 1.5 in.  
Device: L 37.5 in.  
Weight: 2.7 oz



online video available

< scan to view

### C-A-T TOURNIQUET

COLOR	STYLE	ITEM #
Black	Tactical	30-0001
Orange	Rescue	30-0023
Blue	Trainer	30-0033

## COMBAT GAUZE® LE

- Sterile, non-woven Z-fold gauze impregnated with kaolin for temporary external control of bleeding
- Features an embedded x-ray detectable strip



### COMBAT GAUZE LE

ITEM #  
30-0109

### SPECS

Gauze: L 4 yd x W 3 in.  
Pkg: H 5 in. x W 5.5 in. x D 0.625 in.  
Weight: 0.8 oz

## S-ROLLED Gauze

- Unique dispenser, allows controlled linear feed application
- Premium 100% cotton



### SPECIFICATIONS

Gauze: W 4.5 in. x L 4.1 yd  
Packaged: L 3 in. x W 2.5 in. x D 1 in.  
Weight: 1.2 oz

### S-ROLLED

ITEM #  
30-0003



## COMPRESSED Gauze

- Premium rolled 100% cotton gauze for bandaging and wrapping injuries
- Low-cube, ultra-compact sterile

### COMPRESSED

STYLE	ITEM #
Responder	30-0182
Standard	30-0052

### SPECIFICATIONS

Gauze: W 4.5 in. x L 4.1 yd  
Packaged: L 3 in. x W 2 in. x D 1 in.  
Weight: 1.2 oz



## WOUND PACKING Gauze

- Ultra compressed, Z-fold gauze indicated for wound packing
- Soft, white 1-ply gauze with an embedded x-ray detectable line

### WOUND PACKING

ITEM #  
30-0054

### SPECIFICATIONS

Gauze: W 4.5 in. x L 4.1 yd  
Packaged: L 3 in. x W 2 in. x D 1 in.  
Weight: 1.2 oz

## RESPONDER ETD™ Emergency Trauma Dressing

Resilient elastic wrap with sterile non-adherent pad

- Sterile pressure bandages that deliver extremely effective pressure and excellent absorption
- Suitable for applying immediate direct pressure
- Easy-to-use hook & loop securing device with no confusing pressure bars or hooks
- Vacuum-sealed, low-cube packaging



### RESPONDER EMERGENCY TRAUMA DRESSINGS


STYLE	BANDAGE	PAD	WGT	PACKAGE	ITEM #
4 in.Flat	W 4 in. x L 60 in	4x6	2 oz	L 4.25 in. x W 2.75 in. x D 0.88 in.	30-0171
4 in.Flat(Mini)	W 4 in. x L 42 in.	4x6	2.4 oz	L 5 in. x W 3.5 in. x D 0.75 in.	30-0179
6 in.Rolled	W 6 in. x L 60 in.	6x6	2.9 oz	L 6 in. x W 2.5 in. x D 2.1 in.	30-0168



## HYFIN® VENT COMPACT

### Twin Pack

- Two Chest Seals packaged approximately 25% smaller than our standard HyFin®
- 3-channel pressure relief vents
- Two Chest Seals designed for the treatment of entry/exit or multiple penetrating wounds
- Large red tab designed for single step, peel and apply application and for burping the wound if necessary

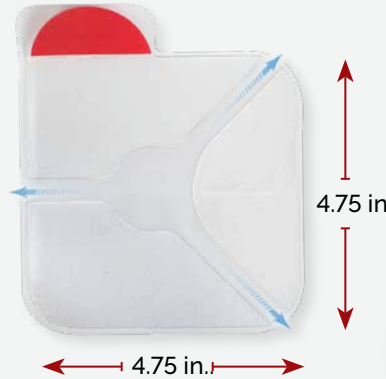

#### HYFIN VENT COMPACT

##### ITEM #

10-0042

#### SPECIFICATIONS

Device (excludes tab): 4.75 in. x 4.75 in.  
Weight: 1.55 oz  
Packaged: H 3.88 in x W 6.6 in x D 0.13 in.



## NASOPHARYNGEAL AIRWAY with lubricant

Packaged with surgical lubricant attached

#### NASO

##### SIZE

##### ITEM #

28F

ZZ-0034

#### SPECIFICATIONS

L 6 in. x W 1.6 in. x D 0.5 in.  
Weight: 0.5 oz

## NAR SURVIVAL BLANKET

Aluminized, reflective PET to help retain up to 90% of your body heat

#### BLANKET

##### ITEM #


ZZ-0038

#### SPECS

Unfolded:  
L 84 in. x W 52 in.  
Weight: 1.7 oz



## Responder NITRILE GLOVES

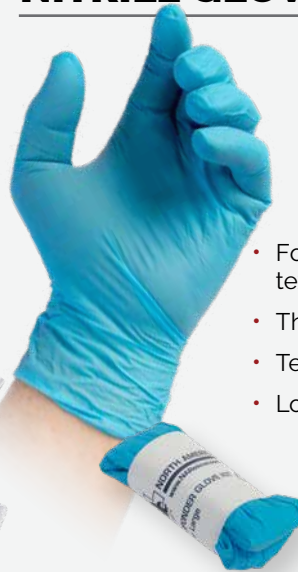

#### SPECIFICATIONS

Kit (25 pair):  
L 4.75 in. x W 6 in. x D 4 in.  
Weight: 1 lb 1 oz



- Form-fitting, ambidextrous design with textured fingertips
- Thick palm ensures durability
- Textured for improved wet and dry gripping
- Longer and thicker than standard nitrile gloves

#### RESPONDER NITRILE GLOVES

SIZE	PKG   QTY	ITEM #
S	Kit   25 pair	ZZ-1059
M	Kit   25 pair	ZZ-1060
L	Kit   25 pair	ZZ-1043
XL	Kit   25 pair	ZZ-1061



## Patient Evacuation | Litters & Stretchers

## NAR QUIKLITTER

LITE, Standard, Responder or XL

#### DISPOSABLE QUIKLITTERS

MODEL	LITTER BED	WGT	PAYLOAD	ITEM #
LITE	L 74 in. x W 27 in.	12 oz	750 lb	60-0071
Standard	L 74 in. x W 28 in.	13.3 oz	500 lb	60-0061
Responder	L 74 in. x W 28 in.	13.3 oz	500 lb	60-0080
XL	L 80 in. x W 40 in.	2.2 lb	1000 lb	60-0065

- Multi-use, low cost disposable evacuation platforms for transport & transfer
- Water-resistant lightweight heavy duty construction with 10 large hand-holds
- Super compact LITE designed to fit in a cargo pocket
- The XL features 12 hand-holds and 2 built-in retention straps





**Active Shooter** incidents  
demand **Proactive Gear**

**MCI  
Rapid  
Response  
Bags** with  
Individual Throw Kits

**THROW KIT**  
INDIVIDUAL KIT WITH COMBAT GAUZE  
KIT CONTENTS:  
■ 1x CAT® Tourniquet, Black  
■ 1x Combat Gauze with Hemostatic Dressing  
■ 1x 4-in. Field Emergency Trauma Dressing  
■ 1x Hybrid Vent Chest Seal, Composite Poly Foam  
■ 1x Cold Blanket/Insulator (Gloves, Large)  
P/N: 80-1105

**C-A-T®** Tourniquet



**NORTH AMERICAN RESCUE®**

[www.NARescue.com](http://www.NARescue.com) • 888.689.6277

**PROTECT** yourself  
& **SAVE** others...

See more lifesaving solutions at  
**[www.NARescue.com](http://www.NARescue.com)**

REV062320