

# Medic Trauma Sheet™



## NAR TACTICAL NYLON SPECIFICATIONS

### FABRIC

We use only the highest quality and most durable fabric available meeting military specification.

Mil-C-43734 C-4

#### **500D CORDURA®**

Solution-Dyed Nylon with Near Infrared Signature Reduction Cordura® nylon is ten (10) times more abrasion resistant than ballistic nylon and the best choice for military load bearing equipment.

A-A 55301 T-3

#### **PALS/MOLLE Webbing**

Solution-Dyed Nylon with Near Infrared Signature Reduction Webbing is an integral part of military load bearing equipment, adding necessary attachment points, suspension and reinforcement. We use only military specification nylon webbing for its strength, durability and chemical resistance.

Mil-PRF-5038 T-3

#### **Binding Tape**

Binding tape protects the cloth edges from fraying. Military load bearing equipment specifications call out a one (1) inch nylon tape that is folded around the cloth edge and sewn in place with a single row of stitching gauged 14 to 5/16 inch from the finished edge.



**NORTH AMERICAN RESCUE®**

www.NARescue.com • 888.689.6277

### THREAD

We use the only suitable thread that has been called out in military load bearing equipment specifications for over 40 years.

V-T-295 T-2 C-A

#### **Thread, Nylon, Size F (92)**

Thread selection is critical to the function and durability of a tactical nylon product. Too light a thread will create weak seams; too heavy a thread will rip the cloth under stress. We use size F (92) nylon thread in all our products for seams and bartacks.

A-A-55126B T-2 C-1

#### **Velcro® Brand Hook & Loop**

Woven nylon Hook 88 & Loop 1000 combination is manufactured and certified to the highest standards

### HARDWARE

We use only hardware meeting rigid military requirements for materials, finish coatings and strength.

N/A

#### **Ghillie-Tex IR Signature Reduction Hardware**

A-A 55634 T-1  
**Reverse Bail  
Slider**

Mil-C-5040

#### **550 Cord Zipper Pull**

NASM 16491 T-1 C-3

#### **Grommet**


NASM 27980

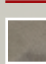
#### **Snap Fastener**


**Rugged, oversized compartments with clear windows for visual organization of your specific contents...**


Designed to provide a versatile, rugged, lightweight equipment access platform for managing casualties during TACEVAC operations. The Trauma Sheet is also an excellent tool for disaster management and trauma rooms. Equipped with rapid access, oversized pockets and large grommet attachment points to facilitate quick set-up and equipment retrieval. Featuring Berry Amendment compliant nylon specs, the Trauma Sheet is constructed of 500D CORDURA® Solution-Dyed Nylon with Near Infrared Signature Reduction, Ghillie-Tex IR signature reduction hardware & silent, self repairing nylon zippers with 550 pull cords.

**Medic Trauma Sheet™**  
COMBAT CASUALTY RESPONSE BAG™

 **ITEM#** 80-0067 (BLACK)  
**NSN#** 6515-01-540-7612

 **ITEM#** 80-0068 (FOLIAGE)  
**NSN#** 6515-01-540-7585

 **ITEM#** 80-0069 (COYOTE)  
**NSN#** 6515-01-540-7603

 **ITEM#** 80-0094 (DUC)  
**NSN#** 6515-01-540-7579

**BAG DIMENSIONS:** Bag Open: H 41.5 in. x  
W 21.5 in. x D 2.5 in.  
Bag Folded: H 20.25 in. x  
W 21.5 in. x D 5 in.  
**WEIGHT:** 4 lb 12 oz



**Medic Trauma Sheet™** ✓  
COMBAT CASUALTY RESPONSE BAG™



**NOTE:**  
Contents are  
not included.



**NORTH AMERICAN RESCUE®**  
[www.NARescue.com](http://www.NARescue.com) 888.689.6277