

HPMK[®] Hypothermia Prevention & Management Kit

REV120222

Top to Bottom Tapered Shape, Maximizes Isothermal Capabilities...

Recommended by CoTCCC for the prevention of hypothermia during casualty care, the patented HPMK[®] with reinforced Heat Reflective Shell is strong, flexible, lightweight & impervious to wind & rain. The HRS[™], the kit's ruggedized outer shell, is constructed of a 4-ply, composite fabric with a protected non-conductive reflective layer that provides excellent thermal insulation. Rapid 360° access to the casualty with minimal exposure to the elements is achieved through the use of 1.5 in. continuous hook & loop closures located along the perimeter of the covering. The HRS[™] also features a built-in hood, fluid absorption pad & top-to-bottom tapered shape to maximize isothermal capabilities. The heat-generating source of the HPMK[®] is a self-heating, oxygen-activated shell liner designed to sustain 10 hours of continuous dry heat (no external power supply required). This low cube kit is vacuum-packaged in a rugged, pliable pouch featuring NAR's Red Tip Technology[®] with tabs that allow easy access and rapid deployment.

HPMK Hypothermia Prevention & Mgt Kit



ITEM#	NSN#
80-0027	6515-01-532-8056



NOTE: Litter is sold separately.




NORTH AMERICAN RESCUE[®]
www.NARescue.com • 888.689.6277

PRODUCTS
WITH  **MISSION[®]**

HPMK[®] HYPOTHERMIA PREVENTION & MANAGEMENT KIT

REV120222

HPMK contents:

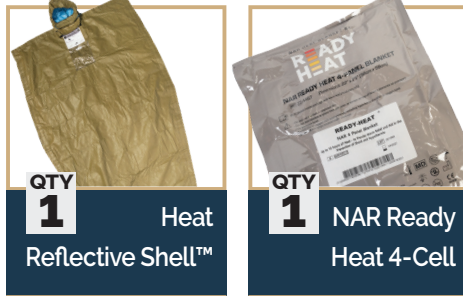

Sustains 10 hours of continuous dry heat

Kit Contents:

- 1 x Heat Reflective Shell (HRS™) with built-in hood
- 1 x Self-heating shell liner
- 1 x Plastic vacuum bag



Features:

- TCCC-recommended hypothermia prevention solution during casualty care
- Sustains 10 hours of continuous dry heat with oxygen-activated, self-heating liner (no external power supply required)
- Constructed of 4-ply composite fabric with protected non-conductive reflective layer that provides excellent thermal insulation
- Tapered shape from top to bottom to maximize isothermal capabilities
- Built-in hood & integrated fluid absorption pad
- Continuous 1.5 in. Velcro® closures for rapid 360° access to the casualty with minimal exposure to the elements
- Reinforced and weight tested to support 250+ lb
- Vacuum-packaged in a durable pouch for a protective, low-cube footprint
- Storage pouch features NAR's Red Tip Technology® with tabs that allow easy access and rapid deployment
- U.S. Patent 7,766,950 B2

Specifications:

- Packaged (HPMK[®]): H 6.75 in. x W 10.5 in. x D 5.5 in.
- Weight (HPMK[®]): 3 lb 8 oz
- Heat Reflective Shell (Open): L 78 in. x W 43 in.

HPMK Hypothermia Prevention & Mgt Kit

ITEM#	NSN#
80-0027	6515-01-532-8056



NORTH AMERICAN RESCUE[®]
www.NARescue.com • 888.689.6277

PRODUCTS
 WITH A
MISSION[®]