

MIGHTYSAT™ Rx OXIMETERS

REV050719



Kit includes MightySat Rx Device, lanyard, 2 AAA batteries and carrying case

MightySat™ Rx Oximeter CE

Size	Item#
Standard	20-0093
Bluetooth LE	20-0094
Bluetooth LE/RRp/PVI	20-0095

Masimo SET® - Proven Technology in Hospitals Worldwide

MightySat™ Rx Fingertip Pulse Oximeter from Masimo features Masimo SET® Measure-through Motion and Low Perfusion™ Pulse Oximetry which overcomes the challenges of low perfusion and motion artifact that limit conventional pulse oximetry. This Signal Extraction Technology® (SET®) uses five parallel signal processing engines that measure through challenging conditions allowing users to accurately measure oxygen saturation and pulse rate even during movement or low blood flow.

Features:

- High Resolution Display - Rotating and adjustable color screen for real-time display of all measurements
- Plethysmographic Waveform - Provides high resolutions display of the patient's pleth waveform
- Touchpad - Provides easy access of settings and device features
- Durable, Lightweight Design - Small and portable for use in clinical and home-care environments
- Flexible Pad - Easily conforms to finger shape for enhanced comfort & secure application to finger to minimize external interference
- Signal IQ (SIQ™) - An indicator of the SpO₂ measurement confidence
- Bluetooth Models (20-0094 & 20-0095) are compatible with free downloadable Masimo Professional Health App

Specifications

- Device: H 2.9 in. x W 1.6 in. x D 1.2 in. (7.4 cm x 4.1 cm x 3.0 cm)
- Display: 3.3 cm
- Weight (with batteries): 2.56 oz (73 g)



NORTH AMERICAN RESCUE®
www.NARescue.com • 888.689.6277

PRODUCTS
 WITH  **MISSION®**

MIGHTYSAT™ Rx OXIMETERS

REV050719

High Resolution Display

- Rotating and adjustable color screen for clear display of all measurements

Plethysmographic Waveform

- Provides depiction of the patient's pleth

Touchpad

- Allows for easy access of settings and device features

Durable, Lightweight Design

- Small and portable for use in clinical environments



	MightySat Rx PN: 20-0093	MightySat Rx with Bluetooth LE PN: 20-0094	MightySat Rx with Bluetooth LE/RRp/PVI PN: 20-0095
Masimo SET® Measure-through Motion and Low Profusion™ Technology	•	•	•
Measurements			
SpO ₂ (Oxygen Saturation)	•	•	•
PR (Pulse Rate)	•	•	•
PI (Perfusion Index)	•	•	•
PVI (Pleth Variability Index)			•
Display			
Pleth Waveform	•	•	•
3.3 cm Color Display	•	•	•
Rotational Screen	•	•	•
Bluetooth LE		•	•
Compatible with Masimo Professional Health App*		•	•
Four Year Limited Warranty	•	•	•



Flexible Pad


- Easily conforms to finger shape for enhanced comfort
- Enables secure application to finger

Free Downloadable Masimo Professional Health App*


Display parameter data to a compatible smart device with Bluetooth® LE enabled models. Measurements can also integrate into Apple's Health app.



VIEW measurements on compatible smart device

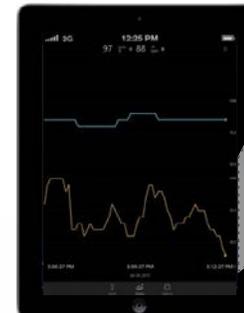

TREND measurements over time, view on graphical display and share the captured data

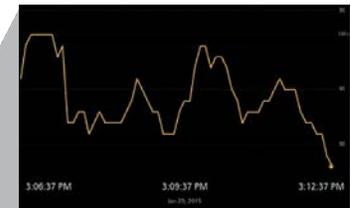

Example of data trend over 6 minutes

*The app is downloadable from the App Store™ for iOS devices or Google Play™ store for select Android devices



NORTH AMERICAN RESCUE®

www.NARescue.com • 888.689.6277

PRODUCTS

WITH

MISSION®