

# EZ-IO® G3 POWER DRIVER NEEDLE SETS (BOX OF 5)

REV071223

For use with the EZ-IO® G3 Power Driver

## Immediate Vascular Access to the central circulation within seconds...

EZ-IO® provides the medical professional immediate vascular access to the central circulation within seconds, delivering medications, intravenous fluids and blood products to adult and pediatric patients. With a specially designed cutting IO needle, manual driver, or small power driver, the EZ-IO® allows the clinician complete control — avoiding the use of force. EZ-IO® provides rapid, smooth entry into the bone's medullary cavity, creating an immediate conduit to the central circulation.

The EZ-IO® Intraosseous Infusion System is the leading choice for intraosseous access across multiple health care settings, and is used by 90 percent of US advanced life support ambulances, over half of US Emergency Departments, and the US Military.

Unlike older IO technology, the EZ-IO®'s patented needle tip has four asymmetrical bevels on the stylet, and two cutting tips on the catheter, providing a smooth and accurate insertion. Black indicator lines are located at 5 mm increments from the hub for assessing insertion depth. With EZ-IO®, insertion depth is completely controlled by the clinician.



15 mm Needle Set\*



45 mm Needle Set\*

**EZ-IO® Products are restricted to Federal Purchases only. Please contact Customer Service for details.**

*\*Caution: Federal (USA) law restricts this device to sale by or on the order of a practitioner licensed to use or order use of the device or Federal (USA) law prohibits dispensing without a prescription.*



25 mm Needle Set\*




**NORTH AMERICAN RESCUE®**  
www.NARescue.com • 888.689.6277

**PRODUCTS**  
WITH A **MISSION®**

*\*(Sold only in Box of 5)*

# EZ-IO® G3 POWER DRIVER NEEDLE SETS (BOX OF 5)

REV071223

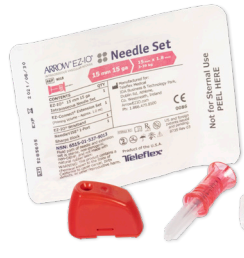

## **EZ-IO® Power Driver Needle Set (Box of 5) Pink Specifications:**

- 15 mm 15 Gauge Needle Set
- Designed for use in patients weighing 3-39 kg (approximately 7 - 86 lb)
- For patients with minimal tissue over insertion sites



## **EZ-IO® Power Driver Needle Set (Box of 5) Blue Specifications:**

- 25 mm 15 Gauge Needle Set
- Designed for any patient weighing more than 39 kg (over 86 lb)
- For patients who have too much tissue over the insertion site for the 15 mm Needle Set to be used



## **EZ-IO® Power Driver Needle Set (Box of 5) Yellow Specifications:**

- 45 mm 15 Gauge Needle Set
- Designed for patients weighing more than 39 kg (over 86 lb)
- For patients with excessive tissue over the targeted insertion site (e.g. edema, large musculature, obesity)
- Ideal for the humerus site in patients over 39kg (over 86 lb)



## **Features:**

- 15 gauge
- Three needle sizes accommodate intraosseous access in patients 3 kilograms and larger
- 304-stainless steel
- Designed for use with any EZ-IO® Power Driver
- NeedleVise® sharps block
- IO wristband (pink for extremity alert)
- EZ-Connect® (IV extension set with a standard Luer connector)
- Black indicator lines are located at 5 mm increments from the hub for assessing insertion depth
- Used by 90 percent of US advanced life support ambulances
- Smooth and accurate insertion

## **Specifications:**

Please see information located next to individual products.

## **Dimensions:**

- (L) 4.5" x (W) 4.5" | Weight (Box): 6.5 oz.

## **EZ-IO® Power Driver Needle Sets**

Color	Size	Item #	NSN#
Pink	15 mm   PD	30-0081	
Blue	25 mm   AD	30-0083	
Yellow	45 mm   LD	30-0085	



STERILE

*The EZ-IO® Intraosseous Infusion System is the leading choice for intraosseous access across multiple health care settings.*



**NORTH AMERICAN RESCUE®**

www.NARescue.com • 888.689.6277

**PRODUCTS**  
WITH A **MISSION®**