

HYFIN® CHEST SEAL DRESSING COMBO PACK

REV070124

The HyFin® Chest Seal Combo Pack from North American Rescue provides one HyFin® Vent Chest Seal and one HyFin® Occlusive Chest Seal packaged in a more compact, rugged, easy-to-open twin foil pouch Combo Pack. Each HyFin® chest seal includes a gauze pad to wipe the wound surface prior to application. The two different seals are each individually packaged and are designed to allow the user to apply a seal to both an entry and exit wound while also giving the option to only apply one as needed and store the other as the situation requires.

The HyFin® Vent Chest Seal's innovative design provides 3-vented channels designed to prevent airflow into the chest cavity during inhalation while allowing air to escape through the vent channels during exhalation. The 3-vent channels are designed to allow blood to escape and also provide a backup fail-safe system, as even if two of the three channels become obstructed, the vent is designed to remain fully operational. The HyFin® Occlusive Chest Seal is an occlusive dressing providing advanced adhesive technology designed to overcome adverse conditions in austere environments where the casualty may be covered in blood, sweat, body hair, or other environmental contaminants when a reliable chest seal is required.

Each chest seal measures 6 in. x 6 in. and has a transparent design that allows for accurate placement directly over the wound, easily conforming to the casualty's chest. A large, easily-visible red elliptical tab, facilitating single step peel-and-adhere application, is incorporated into the HyFin® design. Made in America and 100% latex free, each chest seal is packaged using its signature Red Tip Technology® red tear notches in a rugged, easy-to-open, twin foil pouch to ensure high quality performance. Each side of the HyFin® Chest Seal Combo Kit is clearly identified by name and color to signify which is the HyFin® Vent and which is the HyFin® Occlusive chest seal.



HyFin® Chest Seal Dressing - COMBO

| Item # | NSN# |
|---------|------|
| 10-0073 | ---- |



NORTH AMERICAN RESCUE®
www.NARescue.com • 888.689.6277


PRODUCTS
 WITH A **MISSION®**

Features:

- Versatile Combo Pack designed for the treatment of both entry/exit or multiple penetrating injuries to the chest
- Featuring NAR's patented HyFin® Chest Seals with enhanced adhesive technology designed for a superior seal in the most adverse conditions, including sweaty or hairy casualties
- Each HyFin® measures 6 in. x 6 in. and contains a gauze pad to wipe the wound surface prior to application
- Includes one HyFin® Vent Chest Seal & one HyFin® Occlusive Chest Seal individually packaged in clearly marked and easy-to-open twin foil pouch
- Large, easily-visible red elliptical tabs designed to allow for single step peel-and-apply application and the burping of a wound if necessary

Specifications:

- Device: H 6 in. x W 6 in. (excluding tab)
- Packaged Folded: H 7.5 in. x W 4.5 in. x D 0.25 in.
- Packaged Unfolded: L 7.5 in. x W 9 in. x D 0.13 in.
- Weight: 2.5 oz

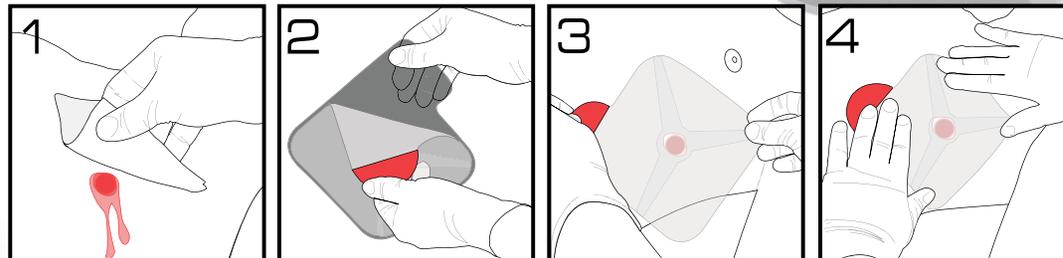

NORTH AMERICAN RESCUE®
www.NARescue.com • 888.689.6277

PRODUCTS
 WITH A **MISSION®**