

KOBUS ROLLUP BULLNOSE LITTER

Learn more about this and other Products with a Mission® at www.NARescue.com



The optimal solution for rescue in confined spaces, heights, depths and horizontal rescues.

The patented RollUP system combines the advantages of a roll-up stretcher with the rigidity of a basket carrier.

The Kobus RollUP Bullnose Litter is designed specifically to optimize the demands required for all aspect of rescue - high angle, confined space, cold weather, marine and fast water applications. Litters are available in 2 different widths: 27.5 and 36.2 inches.

The KOBUS Dura-Slyde (thermo-plastic polymer sliding surface) delivers excellent toughness and strength with a low coefficient of friction. This means no heat generation and outstanding resistance to abrasion and impact during drag maneuvers in a wide range of environmental temperatures and conditions.

KOBUS Litters are designed with 10 stainless steel triangular, permanently affixed lifting / attachment plates with 3 essential grommets per plate for accessories and fixation. Additionally, all Litters feature integrated lifting bridals that make vertical and horizontal hoists a standard capability (no 3rd party straps or adapters needed). Once the patient is secured, the litter becomes rigid. Integrated struts prevent the litter from buckling and prevents pressure being exerted on the patient during the rescue.

KOBUS Litters have been awarded multiple aeromedical certifications for helicopter winching and evacuation operations. All webbing is certified to international standards for lifting, climbing, and flight requirements - including handles and securing straps. All securing straps feature COBRA® steel buckles.

All components have excellent chemical resistance and are easy to clean and decontaminate.

For convenience, the stretcher can be rolled up and stored in an included MOLLE configured backpack and carry case. The carry case also features a removable accessory pouch.

Kohlbrat and Bunz Corporation roll up litters (KOBUS) are **warranted for a period of 10 years from the date manufactured.** Please visit the product page for more information.

FEATURES/SPECIFICATIONS

- Bullnose design is ideal for drag in snow and undergrowth environments, and provides superior head protection in rescue scenarios
- Integrated lifting bridals for vertical and horizontal
- Hoistable. Easy transition between vertical and horizontal hoist
- Fits in confined spaces
- Integrated struts prevent buckling or pressure on the patient
- Material is chemical resistant

DIMENSIONS

27.5 in"

- Product Dimensions: L 96 in. x W 27.5 in.
- Weight: 11 lb.

36.2 in"

- Product Dimensions: L 97.6 in. x W 36.2 in.
- Weight: 18.3 lb.



KOBUS ROLLUP BULLNOSE LITTER

ITEM #	COLOR	WIDTH	NSN
60-0105	Green	27.5 in.	----
60-0107	Green	36.2 in.	----
60-0109	Red	27.5 in.	----
60-0111	Red	36.2 in.	----

* More Kobus products available at www.NARescue.com

Integrated struts on all KOBUS RollUp Litters prevent buckling or pressure on the patient



- 10 stainless steel triangular, permanently affixed, lifting / attachment plates with 3 essential grommets per plate for accessories and fixation are standard
- All webbing is certified to international standards for lifting, climbing and flight requirements including the handles and securing straps (Meets NFPA Standards, European Norm (EN), FDA Class 1 Medical Device)
- Standard Cobra steel buckles on all securing straps
- Multiple aeromedical certifications for helicopter winching / evacuation operations



NORTH AMERICAN RESCUE®
www.NARescue.com • 888.689.6277

KOBUS ROLLUP BULLNOSE LITTER




Bullnose design is ideal for drag in snow and undergrowth environments



Close-up of one of the stainless steel triangular, permanently affixed lifting / attachment plate with 3 grommets



Demonstration of Rollup Enclosed Rescue Stretcher set-up

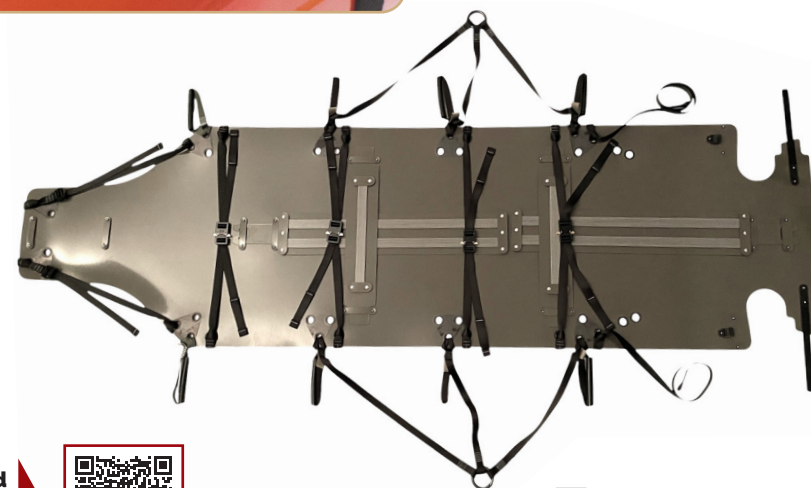

KOBUS ROLLUP BULLNOSE LITTER

ITEM #	COLOR	WIDTH	NSN
60-0105	OD Green	27.5 in.	----
60-0107	OD Green	36.2 in.	----
60-0109	Red	27.5 in.	----
60-0111	Red	36.2 in.	----

* More Kobus products available at www.NARescue.com

Learn more about this and other Products with a Mission® at www.NARescue.com



NORTH AMERICAN RESCUE®
www.NARescue.com • 888.689.6277