

# KOBUS ROLLUP LITTER

Learn more about this and other Products with a Mission® at [www.NARescue.com](http://www.NARescue.com)



**The optimal solution for rescue in confined spaces, heights, depths and horizontal rescues.**

The patented RollUP system combines the **advantages** of a roll-up stretcher with the **rigidness** of a basket carrier.

The Kobus RollUP Litter is designed specifically to optimize the demands required for all aspect of rescue - high angle, confined space, cold weather, marine and fast water applications. Litters are available in 2 different widths: 27.5 and 36.2 inches.

The KOBUS Dura-Slyde (thermo-plastic polymer sliding surface) delivers excellent toughness and strength with a low coefficient of friction. This means no heat generation and outstanding resistance to abrasion and impact during drag maneuvers in a wide range of environmental temperatures and conditions.

KOBUS Litters are designed with 10 stainless steel triangular, permanently affixed lifting / attachment plates with 3 essential grommets per plate for accessories and fixation. Additionally, all Litters feature integrated lifting bridals that make vertical and horizontal hoists a standard capability (no 3rd party straps or adapters needed). Once the patient is secured, the litter becomes rigid. Integrated struts prevent the litter from buckling and prevents pressure being exerted on the patient during the rescue.

KOBUS Litters have been awarded multiple aeromedical certifications for helicopter winching and evacuation operations. All webbing is certified to international standards for lifting, climbing, and flight requirements - including handles and securing straps. All securing straps feature COBRA® steel buckles.

All components have excellent chemical resistance and are easy to clean and decontaminate. For convenience, the stretcher can be rolled up and stored in an included MOLLE configured backpack and carry case. The carry case also features a removable accessory pouch.

Kohlbrat and Bunz Corporation roll up litters (KOBUS) are **warranted for a period of 10 years from the date manufactured**. Please visit the product page for more information.

## FEATURES/SPECIFICATIONS

- Integrated lifting bridals for vertical and horizontal
- Hoistable. Easy transition between vertical and horizontal hoist
- Fits in confined spaces
- Integrated struts prevent buckling or pressure on the patient
- Material is chemical resistant
- Standard Cobra steel buckles on all securing straps
- 10 stainless steel triangular, permanently affixed, lifting / attachment plates with 3 essential grommets per plate for accessories and fixation are standard
- All webbing is certified to international standards for lifting, climbing and flight requirements including the handles and securing straps (Meets NFPA Standards, European Norm (EN), FDA Class 1 Medical Device
- Multiple aeromedical certifications for helicopter winching / evacuation operations



## DIMENSIONS

27.5 in"

- Package Dimensions: H: 10.6 in. x W: 28 in.
- Product Dimensions: L 96 in. x W 27.5 in.
- Weight: 12.1 lb.

36.2 in"

- Package Dimensions: H: 10.6 in. x W: 36.2 in.
- Product Dimensions: L 97.6 in. x W 36.2 in.
- Weight: 17.4 lb.

**Integrated struts on all KOBUS RollUP Litters prevent buckling or pressure on the patient**



## KOBUS ROLLUP LITTER


ITEM #	COLOR	WIDTH	NSN
60-0104	OD Green	27.5 in.	----
60-0106	OD Green	36.2 in.	----
60-0108	Red	27.5 in.	----
60-0110	Red	36.2 in.	----

\* More KOBUS products available at [www.NARescue.com](http://www.NARescue.com)



**NORTH AMERICAN RESCUE®**  
[www.NARescue.com](http://www.NARescue.com) • 888.689.6277

# KOBUS ROLLUP LITTER

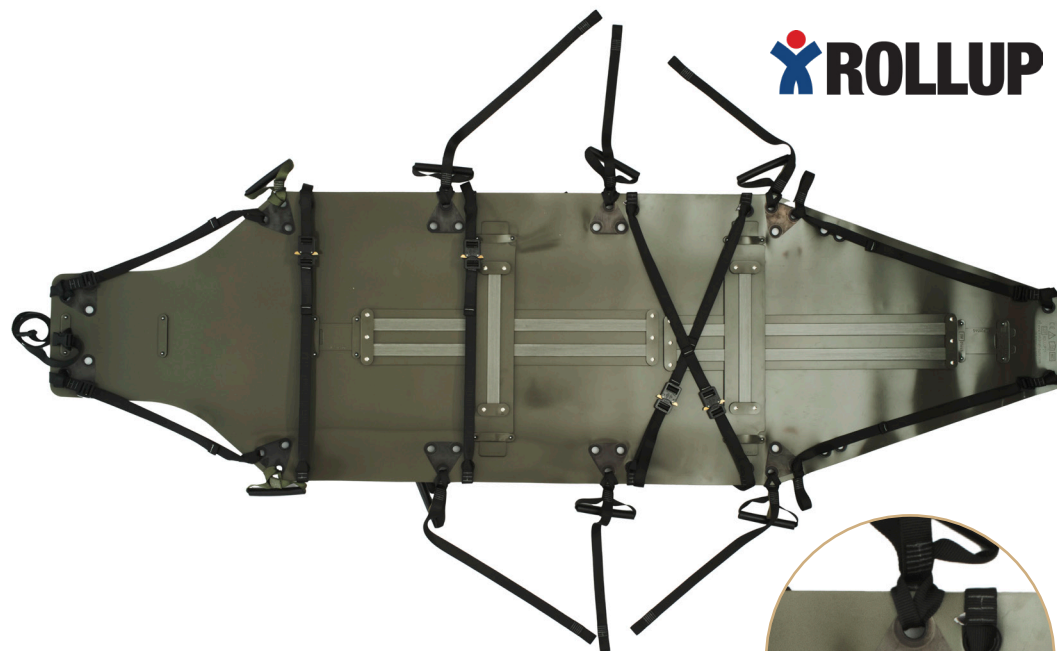

10 steel lifting/attachment plates with 3 grommets



Demonstration of Rollup Rescue Stretcher set-up



## KOBUS ROLLUP RESCUE STRETCHER

ITEM #	COLOR	WIDTH	NSN
60-0104	Green	275 in.	----
60-0106	Green	36.2 in.	----
60-0108	Red	275 in.	----
60-0110	Red	36.2 in.	----

\* More Kobus products available at [www.NARescue.com](http://www.NARescue.com)

Learn more about this and other Products with a Mission® at [www.NARescue.com](http://www.NARescue.com)

