

HEAT REFLECTIVE SHELL - INSULATED™ (HRS-I™)

Learn more about this and other **Products with a Mission**® at www.NARescue.com



The Heat Reflective Shell - Insulated™ (HRS-I™), a rugged outer shell is constructed of a 5-ply composite fabric with a protected non-conductive reflective layer and NEW insulation layer that provides additional thermal retention. Rapid 360° access to the casualty is achieved through the use of 1.5 in. continuous hook and loop closures located along the perimeter of the covering.


The new insulation layer in both the top and bottom of the HRS-I™ exceeds the insulation value of the Army Module Sleep System Patrol Bag. This insulation is highly compressible but also has excellent rebound after opening. The added insulation is also hydrophobic and will not retain fluids or additional liquid weight while maintaining insulation properties.



EXCEEDS ARMY MODULE SLEEP SYSTEM PATROL BAG INSULATION CLO VALUE



5-PLY COMPOSITE FABRIC WITH CLIMASHIELD® INSULATION

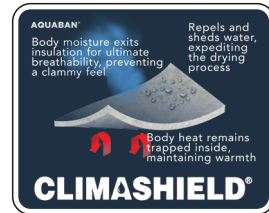

360° RAPID ACCESS HOOK & LOOP CLOSURES

FEATURES/SPECIFICATIONS

- Constructed of 5-ply composite fabric with protected non-conductive reflective layer and NEW insulating layer for improved thermal retention
- The NEW Climashield® insulation layer top and bottom exceed the insulation value of the Army Module Sleep System Patrol Bag in service today. This insulation is highly compressible but also has excellent rebound after opening. The added insulation features Aquaban® technology making it

hydrophobic so that it will not retain fluids or additional liquid weight while maintaining insulation properties.

- Rapid 360° access to the casualty with minimal exposure to the elements is achieved through the use of 1.5 in. continuous hook and loop closures located along the perimeter of the covering.
- Reinforced and weight tested to support 250+ lb

DIMENSIONS

- Packaged: H 11.25 in. x W 9 in. x D 4.25 in.
- Size: L 78 in. x W 43 in.
- Weight: 2.81 lbs

HEAT REFLECTIVE SHELL - INSULATED™

ITEM #	NSN
ZZ-1452	---
ZZ-1452-EU	---



NORTH AMERICAN RESCUE®
www.NARescue.com • 888.689.6277