

# SAM PELVIC SLING

## Ideal compression every time without the guesswork


For rapid, effective stabilization of pelvic fractures. The revolutionary SAM Pelvic Sling is the first and only pelvic belt to contain the patented autostop buckle which is programmed to engage when the correct, safe amount of circumferential pressure is applied. The belt itself is adjustable from 25 in. to 50 in.. MRI safe as long as application on patient is prior to entering the MRI environment.

Learn more about this and other Products with a Mission® at [www.NARescue.com](http://www.NARescue.com)



## FEATURES/SPECIFICATIONS

- The only force-controlled circumferential pelvic belt scientifically & clinically proven to safely & effectively stabilize open-book pelvic ring fractures
- Rapid, simple, three-step application, usually in less than one minute
- Patented Autostop buckle allows a compression force no greater than 33 lb; vital in high-stress environments where over-tightening could potentially be extreme and harmful
- "Click" provides clear feedback to confirm correct application
- Posterior slider for facilitating transfers
- Radiolucent and MRI safe
- Made in America
- Single Use Device
- SEL#: 09ME-04-SPLT

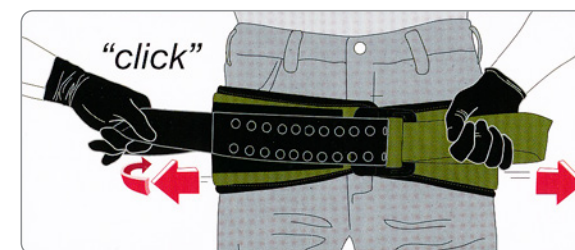

## SAM PELVIC SLING

ITEM #	NSN
50-0004	6515-01-509-6866

## DIMENSIONS

- H 7.25 in. x W 7.75 in. x D 4 in.
- Fits Hip Circumference: from 32 in. to 50 in.
- Weight: 9.8 oz



**NORTH AMERICAN RESCUE®**  
[www.NARescue.com](http://www.NARescue.com) • 888.689.6277