

FIELD CORPSMAN KITS

BASIC & DELUXE

The Basic Field Corpsman Kit from NAR® is a lightweight, durable and portable diagnostic kit provides excellent performance and flexibility in the field. This water resistant, highly durable kit includes a variety of high intensity lights and magnifiers for day or night examinations. Choose from the Basic or Deluxe portable diagnostic kit.

(Deluxe Kit contains all the items from the Basic kit plus a Mini Ophthalmoscope, additional Spare Bulb and extra Batteries).

Fluorescein strips not supplied in kit.



Learn more about this and other Products with a Mission® at www.NARescue.com



FEATURES/SPECIFICATIONS

- Lightweight, portable diagnostic kit ideal for day or night examinations
- 80 lumen LED examination lamp (batteries included)
- Red lens for night use
- Fiber optic light-bender/nasal head
- Spare batteries
- Soft mesh pouch with padded sides and durable zipper
- Fluorescein strips not supplied in kit
- Featuring the Heine Mini 3000 Ophthalmoscope (**Deluxe Only**; batteries included)



Pictured: Deluxe Kit

FIELD CORPSMAN KITS

ITEM #	OPTION	NSN
20-0001	Basic	6515-01-458-6057
20-0002	Deluxe	6515-01-458-6178

* Find individual components at www.NARescue.com.

DIMENSIONS

Deluxe:

L 8 in. x W 4.25 in. x D 2.25 in.

Weight: 11.6 oz

Basic:

• L 8 in. x W 4.25 in. x D 2.25 in.

• Weight: 6.2 oz



NORTH AMERICAN RESCUE®
www.NARescue.com • 888.689.6277


FIELD CORPSMAN KITS

BASIC & DELUXE

all KITS INCLUDE:



+plus, by Kit Type:

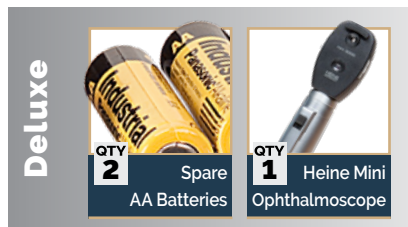

Learn more about this and other Products with a Mission® at www.NARescue.com



BASIC KIT CONTENTS

Kit Contents:

- General examination light: 80 lumen LED lamp to provide a focused beam of bright white light for superior illumination. Lightweight, water resistant and guaranteed unbreakable (batteries included)
- Red lens: Fits on the exam light for night use
- Fiber optic light bender/nasal head: To effectively direct light into the nasal or oral cavity. Can be used to transilluminate the sinuses
- 2 spare AAA batteries
- Soft case with mesh pouches: Designed to secure the exam lights and accessories. Developed with padded sides and durable zipper. Two additional mesh pouches are positioned on the inside compartment for more parts or supplies
- Color instruction booklet: Contains operating instructions and part numbers

Ear Care:

- Otoscopic head: Provides bright illumination and 3X magnification for excellent visualization. Adjustable head allows examiners to adjust for personal comfort
- Specula adaptor: For use with the otoscopic head
- Disposable adult specula: For use with the otoscopic head. Otolight accepts standard disposable specula
- Lighted stainless curette: Used to remove wax from the ear or to assist in debris removal from small wounds

Eye Care:

- Cobalt blue lens: Provides UV illumination to examine the eyes for corneal abrasions or small foreign bodies. Use with or without fluorescent eye dyes
- High resolution eye loop: Soft nylon loop fits on the cobalt blue lens. Designed to gently sweep foreign bodies from the surface of the eye. The lighted tip illuminates work area
- Lighted eye magnet: For the removal of metal foreign bodies in the eye

DELUXE KIT CONTENTS

- All Basic Contents **PLUS:** Heine Mini 3000 Ophthalmoscope (batteries included)



*Pictured:
Basic Kit*

FIELD CORPSMAN KITS

ITEM #	OPTION	NSN
20-0001	Basic	6515-01-458-6057
20-0002	Deluxe	6515-01-458-6178

* Find individual components at www.NARescue.com.



NORTH AMERICAN RESCUE®
www.NARescue.com • 888.689.6277