

KOBUS ROLLUP BULLNOSE LITTER

23.6 IN (60 CM) - GREEN

The Optimal Solution For Rescue In Confined Spaces, Heights, Depths And Horizontal Rescues.

*The patented RollUP system combines the **advantages** of a roll-up stretcher with the **rigidness** of a basket carrier.*

The Kobus RollUP Bullnose Litter is designed specifically to optimize the demands required for all aspect of rescue - high angle, confined space, cold weather, marine and fast water applications. Litters are available in 2 different widths: 27.5 and 36.2 inches.

The KOBUS Dura-Slyde (thermo-plastic polymer sliding surface) delivers excellent toughness and strength with a low coefficient of friction. This means no heat generation and outstanding resistance to abrasion and impact during drag maneuvers in a wide range of environmental temperatures and conditions.

KOBUS Litters are designed with 10 stainless steel triangular, permanently affixed lifting / attachment plates with 3 essential grommets per plate for accessories and fixation. Additionally, all Litters feature integrated lifting bridals that make vertical and horizontal hoists a standard capability (no 3rd party straps or adapters needed). Once the patient is secured, the litter becomes rigid. Integrated struts prevent the litter from buckling and prevents pressure being exerted on the patient during the rescue.

KOBUS Litters have been awarded multiple aeromedical certifications for helicopter winching and evacuation operations. All webbing is certified to international standards for lifting, climbing, and flight requirements - including handles and securing straps. All securing straps feature COBRA® steel buckles.

All components have excellent chemical resistance and are easy to clean and decontaminate.

For convenience, the stretcher can be rolled up and stored in an included MOLLE configured backpack and carry case. The carry case also features a removable accessory pouch.

FEATURES/SPECIFICATIONS

- Bullnose design is ideal for drag in snow and undergrowth environments, and provides superior head protection in rescue scenarios
- Integrated lifting bridals for vertical and horizontal
- Hoistable. Easy transition between vertical and horizontal hoist
- Fits in confined spaces
- Integrated struts prevent buckling or pressure on the patient
- Material is chemical resistant
- Standard Cobra steel buckles on all securing straps

DIMENSIONS

23.6 in*

- Product Dimensions: L 96 in. x W 23.6 in.
- Weight: 11 lb..
- Weight (with carrying case): 15.5 lb.
- Weight (without carrying case): 14.3 lb.

Learn more about this and other Products with a Mission® at www.NARescue.com



ROLLUP

Integrated struts on all KOBUS RollUp Litters prevent buckling or pressure on the patient

KOBUS ROLLUP BULLNOSE LITTER

ITEM #	COLOR	WIDTH	NSN
60-0117	Green	23.6 in.	----

* More Kobus products available at www.NARescue.com

- 10 stainless steel triangular, permanently affixed, lifting / attachment plates with 3 essential grommets per plate for accessories and fixation are standard
- All webbing is certified to international standards for lifting, climbing and flight requirements including the handles and securing straps (Meets NFPA Standards, European Norm (EN), FDA Class 1 Medical Device)
- Multiple aeromedical certifications for helicopter winching / evacuation operations



Inflatable Pillow Sold Separately



NORTH AMERICAN RESCUE®
www.NARescue.com • 888.689.6277

KOBUS ROLLUP BULLNOSE LITTER

23.6 IN (60 CM) - GREEN




Demonstration of Bullnose Rollup Rescue Stretcher Set-up

KOBUS ROLLUP BULLNOSE LITTER

ITEM #	COLOR	WIDTH	NSN
60-0117	Green	23.6 in.	----

* More Kobus products available at www.NARescue.com

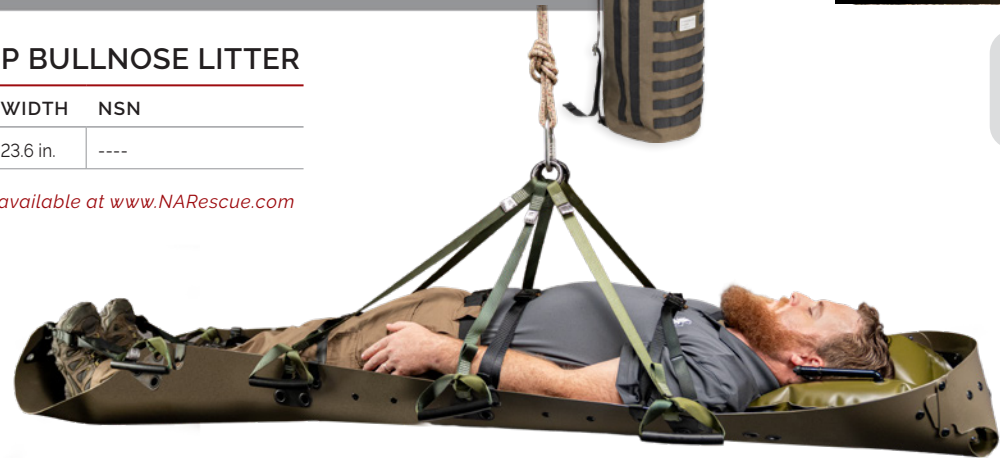

Learn more about this and
other Products with a Mission®
at www.NARescue.com



Close-up of one of the stainless
steel triangular/triangular
permanently affixed lifting/attachment
plates with 3 grommets

Kohlbrat and Bunz Corporation roll up litters (KOBUS) are
warranted for a period of 10 years from the date manufactured.
Please visit the product page for more information.

Inflatable Pillow
Sold Separately



NORTH AMERICAN RESCUE®
www.NARescue.com • 888.689.6277