

LifeFlow RAPID INFUSER *from 410 Medical*

REV121720

Deliver 500ml of critical fluids in less than 2 minutes

Born out of the barriers that users often face delivering fluids quickly and efficiently to critically ill patients, LifeFlow is a hand-operated rapid infusion device that allows health care providers to deliver fluids quickly and efficiently, improving outcomes for critically ill patients. Timely fluid resuscitation is key to the effective treatment of severe sepsis and shock, but is an often-neglected component of care due difficulties with vascular access and the slow and inefficient fluid delivery techniques commonly used.

LifeFlow comes in a simple kit that's easy to set up and operate. The one-handed operation lets you work on more than one critical task at a time. LifeFlow makes emergency fluid resuscitation faster, simpler, and more controlled.

Features:

- **Fast:** Deliver 500ml in less than 2 minutes
- **Easy to Use:** Set up and ready for use in less than 3 minutes by a new user
- **Convenient:** Light and portable for use in almost any environment. Disposable and never needs cleaning
- **Safe:** Designed to protect the syringe from bacterial contamination

LifeFlow RAPID INFUSER COMPONENTS

ITEM #	COMPONENTS
30-0286	LifeFlow Fluid Infuser
30-0287	LifeFlow PLUS Blood & Fluid Infuser
30-0288	LifeFlow PLUS Replacement Tubing Set



LifeFlow®
VOLUME RESUSCITATION. WHEN MINUTES MATTER.



NORTH AMERICAN RESCUE®

www.NARescue.com • 888.689.6277

PRODUCTS
— WITH —
MISSION®

LifeFlow RAPID INFUSER *from 410 Medical*

REV121720

SPEED

Average time to set up and prime LifeFlow is less than 2 minutes, and LifeFlow can deliver 500 ml of fluid in less than 2.5 minutes through a 20G IV catheter.

CONTROL


LifeFlow delivers fluid in 10mL increments with each handle compression. This syringe size allows the user to overcome the resistance of IV catheters. The syringe automatically refills with handle release, allowing the user full control of volume of fluid delivered.

EFFICIENCY

LifeFlow can be operated with one hand, leaving the operator's other hand free to address other clinical issues, monitor the infusion site, and assess patient response. This is particularly useful in pre-hospital, emergency department, and ICU settings where speed and efficiency are important and personnel resources may be limited.

LifeFlow®

VOLUME RESUSCITATION. WHEN MINUTES MATTER.

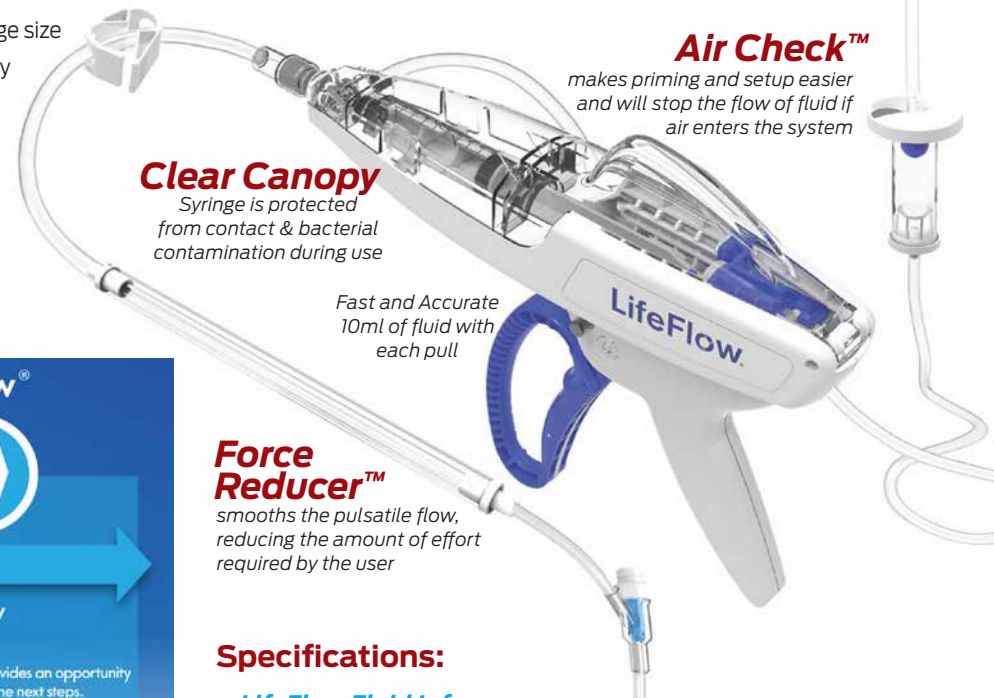

Air Check™

*makes priming and setup easier
and will stop the flow of fluid if
air enters the system*

Clear Canopy

*Syringe is protected
from contact & bacterial
contamination during use*


*Fast and Accurate
10ml of fluid with
each pull*

Force Reducer™

*smooths the pulsatile flow,
reducing the amount of effort
required by the user*

Specifications:

- **LifeFlow Fluid Infuser**
Weight: 0.95 lb
- **LifeFlow PLUS Blood & Fluid Infuser**
Weight: 0.95 lb
- **LifeFlow PLUS Replacement Tubing Set**
Weight: 0.34 lb

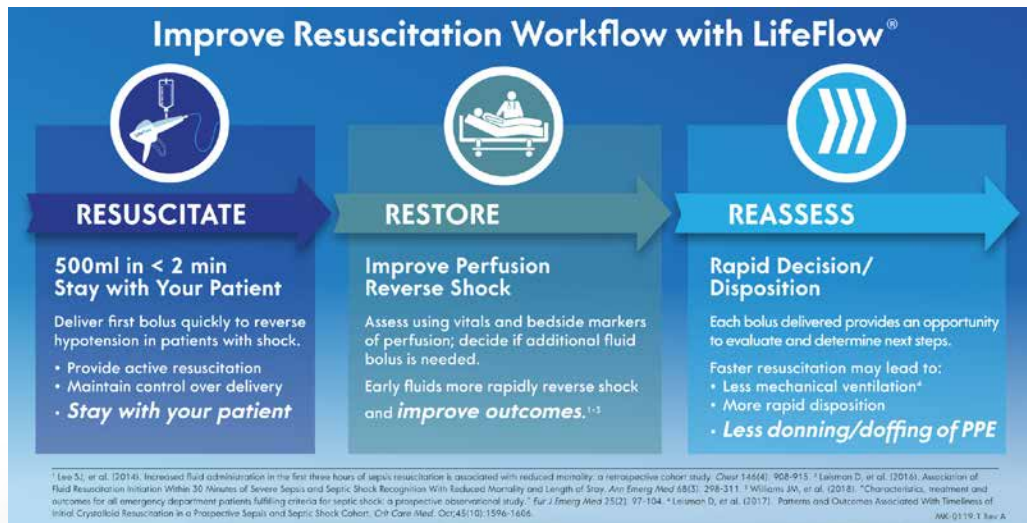

NORTH AMERICAN RESCUE®

www.NARescue.com • 888.689.6277

PRODUCTS
— WITH —
MISSION®