

Ballistic Protection



BALLISTIC PROTECTION | BALLISTIC VEST SYSTEMS | MEDICAL SUPPLIES

OVERVIEW | Body Armor

In recent times there has been a dramatic increase in violent crime and an alarming rise of attacks directed at first responders. Now more than ever, first responders need to consider their personal safety and adopt body armor as PPE. Armor has typically been associated with tactical medics in high threat environments but first responders are realizing the critical need for protection on a daily basis.


Fire and EMS have traditionally been limited to adapting body armor originally designed for the Military & Law Enforcement, North American Rescue is now confronting this challenge by designing ballistic solutions specifically for the unique needs of the Fire & EMS community.

Soft & Hard | 2 Main Types of Body Armor

SOFT ARMOR offers critical protection against handgun threats along with the capability for extended daily wear. NAR provides an optimal over-the-uniform soft armor vest that features flexible layers of ballistic resistant materials constructed from synthetic Aramid fibers with a very high strength to weight ratio. These super strong fibers are woven into multiple layers that are stacked to create an armor panel. When a handgun bullet strikes the armor panel, it gets caught in a web of high performance fibers that stretch and disperse the impact energy causing the bullet to deform or "mushroom". This dispersion of energy not only prevents bullet perforation, but also helps minimize possible injury to underlying tissue.

HARD ARMOR is designed to provide a higher degree of protection from rifle threats. Constructed from a variety of hard ballistic-resistant materials such as ceramic, compressed UHMWPE (ultra-high-molecular-weight-polyethylene), metallic plates, or a composite of multiple materials, each of these rigid materials offer their own unique properties. The hard armor plates offered by NAR are stand-alone (no soft armor required) and multi-hit rated for armor piercing threats such as 5.56 mm "Green Tip" (M855). ▶

How Soft Armor disperses energy from a handgun projectile...

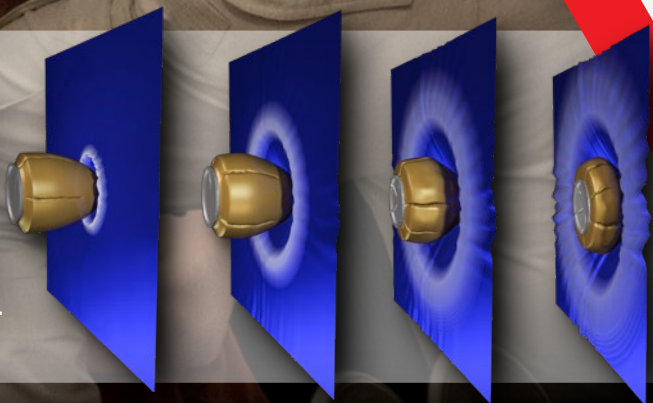

Handguns are overwhelmingly used in most violent crime while Rifles are used 58% of the time in active shooter situations.

Level IIIA soft armor is the foundation of ballistic PPE providing the most coverage against handguns. Rifle plates should be considered for use in any high threat environment.



NIJ 0101.06 Ballistic Standards

How to Choose the right **ARMOR** for YOU...

There are several key features to consider in choosing the right body armor: **Protection Level | Weight | Comfort | Mobility**, and **Cost**. It is important to choose an armor system that will provide protection but also enhance your ability to perform your job. Particularly for Fire & EMS personnel, it is critical that the system facilitates your ability to provide life saving skills and equipment.

North American Rescue® offers the optimal ballistic solution designed as a system specifically for Fire & EMS personnel and the duties they must perform.

The **National Institute of Justice (NIJ)** is responsible for establishing performance standards and testing certification for body armor. The NIJ conducts two types of ballistic testing: Perforation/BackFace Signature (P-BFS) and Ballistic Limit (V50). No complete or partial perforation of the armor panel is allowed during testing. In addition, no more than 44mm (1.73 inches) of indentation is allowed into the ballistic clay backing utilized during shoot testing. The back face signature translates to the upper tolerable limits of blunt force injury to the user. The ballistic limit (V50) is the velocity at which the projectile is expected to perforate the armor 50% of the time. It is important to ensure that any armor package has been certified by an NIJ approved laboratory.

Special Threats

5.56mm "Green Tip" (M855)

was designed for military use and has a steel core. It is often the cheapest and most widely available AR-15 ammo on the market today. Level III plates should be NIJ special threat tested to stop the green tip.



Level II:

Tested to stop 9mm and .357 Magnum ammunition

Level IIIA:

Tested to stop .357 SIG and .44 Magnum ammunition

Level III:

Tested to stop 7.62mm FMJ rifle ammunition

Level IV:

Tested to stop .30 cal steel core armor piercing rifle ammunition

◀ **NIJ Threat Projectiles & Related Armor Levels**



Ballistic Armor

NAR® RESPONDER BALLISTIC PPE VEST

WITH MOLLE & Armored Cummerbund

- Designed specifically for Fire and EMS personnel providing Level IIIA soft armor protection & MOLLE-style PALS attachment capability
- Lightweight and comfortable for a better fit that allows full range of motion
- NIJ 0101.06 Level IIIA ballistic panels designed for maximum coverage without bulk
- Fitted external plate pocket to increase comfort
- Includes Armored Cummerbund with hook and loop adjustment via front/rear kangaroo panel to reduce armor overlap. Small cummerbund option sold separately
- Optional armored lower abdomen extension provides additional coverage with MOLLE-style PALS webbing for essential gear
- Rapidly adjustable shoulder straps
- 3" x 10" Hook & Loop areas on the front & back to accommodate ID Patches (sold separately)
Contact customer service for details.
- PFAS Free

SOFT ARMOR BALLISTIC SPECS

Protection Type/Threat Level:
NIJ 0101.06 IIIA Special Threat Tested
Construction: Woven Twaron / UD Twaron
Thickness: 0.25 in. | AD = 118 psf
NIJ certified Model#: SLSIIIA

NIJ 0101.06 Testing Certification:

V50 - 357sig 125gr TMJ (new) | 1867 ft/s
P-BFS Avg. - 357sig 125gr TMJ | 27.3 mm
V50 - 44mag 240gr JHP (new) | 1661 ft/s
P-BFS Avg. - 44mag 240gr JHP | 36 mm

Special Threats:

V50 - 2gr RCC | 3060 ft/s
V50 - 4gr RCC | 2642 ft/s
V50 - 16gr RCC | 2244 ft/s
V50 - 17gr FSP (.22cal) | 2017 ft/s
V50 - 64gr RCC | 1866 ft/s
V50 - 9mm 127gr SXT (new) | 1754 ft/s
P-BFS Avg. - 9mm 127gr SXT | 24.6 mm
P-BFS Avg. - 40cal 165gr S+W Sp GDHP | 27.2 mm
P-BFS Avg. - 5.7x28 27gr SS195 | 17 mm
P-BFS Avg. - 5.7x28 27gr SS197 | 19.3 mm

SPECIFICATIONS

Vest: L 15 in. x W 18 in. x D 0.5 in..
Waist Size expandable up to 62 in.
Total Weight: 9.1 lb



▲
Vest/Cummerbund options include Reflective Tape Trimmed with Yellow or Orange



ID Patches Sold Separately



NAR® RESPONDER BALLISTIC VEST WITH MOLLE & Cummerbund

COLOR TRIM	ITEM #
Black	70-0331
Navy	70-0333
Red	70-0334
Coyote	70-0332
Multicam	70-0379
Blk Org Reflect	70-0335
Blk Yel Reflect	70-0336
Coy Org Reflect	70-0337
Coy Yel Reflect	70-0338
Nav Yel Reflect	70-0339
Red Org Reflect	70-0340
Red Yel Reflect	70-0341

For New York customers, purchase of this product is limited to eligible professions. Please contact Customer Service with any questions.

NAR® BALLISTIC PPE CARRIER

- Store ballistic vest, plates and helmet
- Highly mobile and easily accessible
- Durable nylon construction
- Nylon webbed carrying strap and handles
- Clear window for holding identification card

SPECIFICATIONS

L 23 in. x W 23 in. x D 8 in.
Weight: 1.43 lbs.



NAR® BALLISTIC PPE CARRIER

ITEM #
70-0372

Ballistic Armor

NAR® RESPONDER BALLISTIC PPE VEST

WITH Armored Cummerbund

- Cost effective option of our NIJ 0101.06 Level IIIA ballistic vest designed for maximum coverage without bulk
- Lightweight, comfortable, better fit allows full range of motion while also providing superb coverage
- Fitted external plate pocket to increase comfort
- Includes Armored Cummerbund with hook and loop adjustment via front/rear kangaroo panel to reduce armor overlap
- PFAS Free
- Rapidly adjustable shoulder straps
- 3" x 10" Hook & Loop areas on the front & back to accommodate ID Patches (sold separately)
- Small cummerbund option is sold separately
- Optional armored lower abdomen extension (sold separately) provides additional coverage with MOLLE-style PALS webbing for essential gear

SPECIFICATIONS

Vest: L 15 in. x W 18 in. x D 0.5 in.
Waist Size expandable up to 62 in..
Total Weight: 8.96 lb

NAR® RESPONDER BALLISTIC VEST WITH Cummerbund

COLOR TRIM	ITEM #	Blk Yel Reflect	70-0364
Black	70-0359	Coy Org Reflect	70-0365
Navy	70-0361	Coy Yel Reflect	70-0366
Red	70-0362	Nav Yel Reflect	70-0367
Coyote	70-0360	Red Org Reflect	70-0368
Blk Org Reflect	70-0363	Red Yel Reflect	70-0369

For New York customers, purchase of this product is limited to eligible professions. Please contact Customer Service with any questions.



▲ Abdominal Extender can be used on both MOLLE and Slick versions (Vest sold separately)



NAR® Responder Armor ABDOMINAL EXTENDER

- Same superior ballistic protection as our NIJ level IIIA soft armor
- Provides additional ballistic protection and platform to mount life saving gear
- Attaches via hook and loop and allows for full range of motion
- PFAS Free

SPECIFICATIONS

H 4.5 in. x W 10 in. x D 0.5 in.
Weight: 0.38 lb

ABDOMINAL EXTENDER

COLOR	ITEM #
Black	70-0342
Navy	70-0344
Red	70-0345
Coyote	70-0343

For New York customers, purchase of this product is limited to eligible professions. Please contact Customer Service with any questions.

NAR® Responder **SMALL CUMMERBUND**

SPECIFICATIONS

L 22 in. x W 6 in. x D 0.5 in.
Maximum Waist Size: 52 in.
Weight: 1.71 lb

SMALL CUMMERBUND

COLOR TRIM	ITEM #
Black	70-0346
Navy	70-0348
Red	70-0349
Coyote	70-0347
Blk Org Reflect	70-0350
Blk Yel Reflect	70-0351
Coy Org Reflect	70-0352
Coy Yel Reflect	70-0353
Nav Yel Reflect	70-0354
Red Org Reflect	70-0355
Red Yel Reflect	70-0356

CHEST POUCH MED KIT

- Designed to add medical equipment to RTF Ballistic Vests with MOLLE-style attachment
- Configurable based on individual skill set and mission requirements
- Rapid access & organization enhanced by elastic storage loops that keep items in place until needed
- Sand and water resistant zippers with silent pulls
- PALS/MOLLE-style attachment system
- Large, ambidextrous Trauma Shear configuration



SPECIFICATIONS

H 7.75 in. x W 9.75 in. x D 4.5 in
Weight: 2 lb 10.3 oz


CHEST POUCH MED KIT



QTY 2
C-A-T®
Tourniquet



QTY 2
ETD® 6 in.
Responder

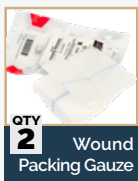

QTY 2
Wound
Packing Gauze



QTY 2
Pre-lubed
NAR® Naso, 28F



QTY 1
HyFin®
Vent Compact, Twin



QTY 1
Responder
Trauma Shears, Lg



QTY 1
Scissor
Leash™



PAIR 3
Responder
Nitrile Gloves Lg



QTY 4
Combat
Casualty Card



QTY 2
Permanent
Marker, Large

CHEST POUCH MEDICAL KIT

COLOR	ITEM #
Black	80-0407
Coyote	80-1109
Red	80-0406

Custom ID PATCHES




Contact customer service to specify your preference.

Options include colors, titles and lettering (solid or reflective)

Ballistic Armor System


NAR® RESPONDER BALLISTIC VEST SYSTEM WITH Medical Kit

COLOR TRIM	ITEM #
Black	80-1213
Navy	80-1215
Red	80-1216
Coyote	80-1214
Multicam	80-1534
Blk Org Reflect	80-1217
Blk Yel Reflect	80-1218
Coy Org Reflect	80-1219
Coy Yel Reflect	80-1220
Nav Yel Reflect	80-1221
Red Org Reflect	80-1222
Red Yel Reflect	80-1223

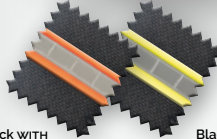

Black WITH
Orange | Reflective


Black WITH
Yellow | Reflective

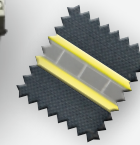

Navy WITH
Yellow | Reflective



Custom ID
Patches
for display
purposes only,
Each sold
separately

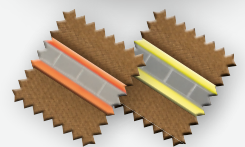

Coyote WITH
Orange | Reflective

Coyote WITH
Yellow | Reflective

SOFT ARMOR BALLISTIC SPECS


Protection Type/Threat Level:
NIJ 0101.06 IIIA Special Threat Tested
Construction: Woven Twaron / UD Twaron
Thickness: 0.25 in. | AD = 118 psf
NIJ certified Model#: SLSIIIA

NIJ 0101.06 Testing Certification:

V50 - 357sig 125gr TMJ (new) | 1867 ft/s
P-BFS Avg. - 357sig 125gr TMJ | 27.3 mm
V50 - 44mag 240gr JHP (new) | 1661 ft/s
P-BFS Avg. - 44mag 240gr JHP | 36 mm

Special Threats:

V50 - 2gr RCC | 3060 ft/s
V50 - 4gr RCC | 2642 ft/s
V50 - 16gr RCC | 2244 ft/s
V50 - 17gr FSP (22cal) | 2017 ft/s
V50 - 64gr RCC | 1866 ft/s
V50 - 9mm 127gr SXT (new) | 1754 ft/s
P-BFS Avg. - 9mm 127gr SXT | 24.6 mm
P-BFS Avg. - 40cal 165gr S+W Sp GDHP | 27.2 mm
P-BFS Avg. - 5.7x28 27gr SS195 | 17 mm
P-BFS Avg. - 5.7x28 27gr SS197 | 19.3 mm

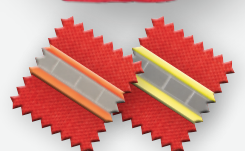

Red WITH
Orange | Reflective

Red WITH
Yellow | Reflective

Vest/Cummerbund options include Trim in Yellow or Orange with Reflective Tape

NAR® RESPONDER ^{WITH Medical Kit} BALLISTIC VEST SYSTEM

The optimal solution for Fire, EMS and RTF personnel in active shooter and mass casualty incidents, with an emphasis on integrating Tactical Emergency Casualty Care...

- NIJ 0101.06 Level IIIA Soft Armor Ballistic Protection equipped with medical supplies that address the leading causes of preventable death in the tactical environment
- Fitted external plate pockets designed for added comfort to accommodate additional hard armor plates (each sold separately)
- Includes armored cummerbund for increased side protection with easy hook and loop attachment for rapid user adjustment expandable waist size up to 62 in.
- 3" x 10" Hook & Loop areas on the front & back to accommodate ID Patches (sold separately)
Contact customer service for details.

PFAS Free

- NAR® Responder Ballistic Vest System includes:
 - 1 x NAR® Responder Ballistic Vest with Armored Cummerbund
 - 1 x Chest Pouch Medical Kit with Luminous Med ID Patch
 - 2 x Nylon Tourniquet Holder with C-A-T®



Custom ID Patches for display purposes only. Each sold separately

SPECIFICATIONS

Vest, Cummerbund & Medical Kit:
L 15 in. x W 18 in.
x D 5 in.

Total Weight:
12.4 lb

Waist Size
expandable
up to 62 in.

Medical Kit
(Stand alone):
H 7.75 in. x W 9.75 in.
x D 4.5 in
Weight: 2 lb 10.3 oz




Chest Pouch Medical Kits are compatible with NAR® Quick Detach Panels

on CUMMERBUND:

CHEST POUCH MEDICAL KIT contains:



▲ NAR®
Quick
Detach
Panel
Compatible
(9 in. x 6.5 in.)

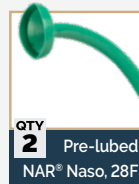

◀ Fully
equipped
Medical Chest
Pouch & TQ's
with Holders
included

NIJ Level III+ STEEL PLATE

- NIJ 0101.06 Level III Certified with Special Threats
- 13/64 or 5mm thickness
- Polyurea coating
- 10 year warranty
- Maximum durability
- Multi-hit Rifle Defense Capability
- Single Curve
- Made in the USA
- **PFAS Free** - made without polyfluoroalkyl substances

SPECIFICATIONS

W 10 in. x H 12 in. x D 0.2 in.
Weight: 6.9 lb

Special Threat Performance:

- 5.56mm x 45mm M193 stops up to 3,000 FPS
- 5.56mm x 45mm M855 / SS109 steel penetrator (Green Tip)
- 7.62mm x 39mm 123 gr steel core
- 7.62mm x 51mm M80 NATO (.308)
- All types of handguns

NIJ LEVEL III+ STEEL PLATE

ITEM #

70-0357

For New York customers, purchase of this product is limited to eligible professions. Please contact Customer Service with any questions.



NIJ Level III+ UHMWPE PLATE Multi Curve - Shooter Cut



***Tested to stop
M855 5.56 "Green Tip"**
- often the cheapest, most
widely available AR15
ammunition



Special Threat Performance:

- 5.45 x 39mm, 53gr (MSC) 7N6
- 5.45 x 39mm, 55gr 7N10
- 5.56 x 45mm, 55gr M193 Ball
- 5.56 x 45mm, 62gr, M855 Green Tip
- 7.62 x 39mm, 123gr, MSC
- 7.62 x 39mm 121gr BZ API
- 7.62 x 54Rmm 148gr LPS

NIJ LEVEL III+ UHMWPE PLATE

ITEM #

70-0428

For New York customers, purchase of this product is limited to eligible professions. Please contact Customer Service with any questions.

- Made with Ceramic and Durable UHMWPE (Ultra High Molecular Weight Polyethylene) materials which are processed to create a lighter, thinner, high performing plate
- High Performing NIJ 0101.06 - Certified Level III
- User-friendly design fits to the contours of the body for maximum safety and comfort
- Rugged, reliable, and stress fracture resistant
- Fully sealed to protect against liquids, chemicals, and harsh elements
- Multi-hit Rifle Defense Capability
- Shooter's Cut Multi Curve
- Stand Alone
- 10 year warranty
- **PFAS Free** - made without polyfluoroalkyl substances

SPECIFICATIONS

W 10 in. x H 12 in. x D 1 in.
Weight: 4.3 lb

WX Saber Advanced BALLISTIC GLASSES

SPECIFICATIONS

FRAME SPECS:

Size: Bridge - 25 mm
Temple - 125 mm
Width - 44 mm
Fit: Universal
Color: Matte Black

BALLISTIC LENSES:

Colors: Smoke Grey & Clear
Ballistic Standards:
MIL-PRF-32432(GL)
ANSI Z87.1-2010
OSHA 1910.133(b)(1)(1)



- Shatterproof Selenite™ Polycarbonate ballistic lenses that meet the MIL-PRF-32432(GL) ballistic standards: ANSI Z87.1-2010 high velocity and high mass impact standards, as well as US Federal OSHA 1910.133(b)(1)(1) standards
- Clear lenses provide maximum light transmission allowing wearer's to see precise color values. Perfect for hazy or overcast, dusk/dawn conditions and indoor use. Light transmission: 89%
- Smoke grey lenses provide glare reduction without distorting colors. Excellent in bright/glare conditions. Light transmission 14%-20%
- T-Shell™ Lens coating resists scratching
- 100% UVA/UVB protection with distortion free clarity
- Adjustable take flight nosepiece accommodates any face/ nose bridge
- Features thin temples compatible with hearing protection
- Includes the PTX UV/Blue/Green Laser Lens Insert to block a variety of lasers covering the light spectrum and provide protection against the largest laser threats

WILEY X SABER ADVANCED BALLISTIC GLASSES

OPTION	ITEM #
Clear & Smoke w/ PTX Insert	90-0108
Clear & Smoke Lenses	90-0114
Clear Lens	90-0109
PTX Insert only	90-0107

PTX LASER LENS INSERT



SPECIFICATIONS

Laser Protection:
VIOLET 400nm OD4+
BLUE 445nm OD4+
GREEN 532nm OD4+

Visible Light
Transmission: 30.7%

Insert Carrier Size:
Width - 44 mm
Bridge - 25 mm

PTX LASER LENS INSERT

ITEM #
90-0107

- **Targeted Laser Protection**
Blocks a variety of lasers covering the light spectrum (**UV/Blue/Green**) to provide protection against the largest laser threats
- Fits Saber Ballistic Glasses
- Frame/Carrier sold separately

EPIC™ RESPONDER PLUS BALLISTIC HELMET

Team Wendy®

- Substantially lighter than standard ACH helmets
- Offers maximum adjustability and comes configured with the EPIC Air liner system, Cam Fit retention and front cover plate
- Shell is constructed from Ultra High Molecular Weight Polyethylene (UHMWPE) and features no thru-hole construction for retention system
- 10 year warranty on ballistic shell from date of delivery
- Tested to NIJ IIIA (according to NIJ STD 0106.-01/0108.01)
- Team Wendy's EPIC™ Responder line of ballistic helmets features lightweight high-performance material paired with the Team Wendy liner systems for premium comfort.



PERFORMANCE STANDARDS

Ballistic Performance: Tested to NIJ IIIA (according to NIJ STD 0106.-01/0108.01)

Fragment Performance: 17 gr V50 > 780 m/s

Backface Deformation: 9mm BFD < 25.4 mm

Blunt Impact Performance: Meets ACH Gen II Blunt Impact (AR/AP 1401) and AEP 2902 clause 4.111 Method H

SPECIFICATIONS:

Head Circumference:
Sm: 20.5 - 21.25 in.
Md: 21.26 - 22.4 in.
Lg: 22.5 - 23.5 in.
Xl: 23.6 - 25.2 in.

EPIC™ Responder Plus Ballistic Helmet

COLOR	SIZE	ITEM#	COLOR	SIZE	ITEM#
● BLK	SM	70-0416	● RGR	SM	70-0420
● BLK	M	70-0417	● RGR	M	70-0421
● BLK	LG	70-0418	● RGR	LG	70-0422
● BLK	XL	70-0419	● RGR	XL	70-0423

ADVANCED COMBAT HELMET *Full Cut*

- Protection Level IIIA according to NIJ 0106.01
- Tested and verified to meet or exceed ballistic resistance as specified under NIJ Standard 0106.01 by an NIJ accredited laboratory
- Materials Including: Kevlar® Fibers and thermoset resin shell construction
- Head Protection: Front, back, crown, and top of head
- Harness: Dial ratchet option and comes standard with 4 pt Harness
- Padding: D30 TRUST 7 Pad helmet system



BALLISTIC SPECS

NIJ 0101.06 Testing Certification:

357 sig | 125 gr FMJ | 1470 ft/s

44 mag | 240 gr SJHP | 1430 ft/s

357 mag | 158 gr JSP | 1430 ft/s

44 mag | 240 gr SWCGC | 1400 ft/s

9 x 19 mm Parabellum | 124 gr FMJRN | 1400 ft/s

40 S&W | 180 gr FMJ | 1155 ft/s

SPECIFICATIONS:

Dimensions:
H 7.9 in. - 8.3 in.
W 6.3 - 6.5 in.
Cir. 22.7 in. - 24.0 in.
Weight: 3 lbs.

SS 401 Advanced Combat Helmet - Full Cut

COLOR	SIZE	ITEM#	COLOR	SIZE	ITEM#
● BLK	SM	70-0390	● TAN	SM	70-0394
● BLK	M	70-0389	● TAN	M	70-0393
● BLK	LG	70-0388	● TAN	LG	70-0392
● BLK	XL	70-0391	● TAN	XL	70-0395



HASTY HARNESS

Multi-purpose Harness

- Indispensable gear for every webbing equipment toolbox
- Resourceful, multi-purpose harness for unlimited applications & maximum versatility
- Constructed of 22 ft loop of 1 in. Tech Tape (MBS: 4,400 lb) with opposing color thread to aid in visual safety inspection

SPECIFICATIONS

Loop: L 10 ft 11 in. x W 1 in.
Circumference - 21 ft 10 in.
Min. Break Str: Tech Tape (4,400 lb)
Sewn Sling (6,000 lb)

HASTY HARNESS

COLOR	ITEM #
Red	90-0044
Black	90-0001



INDIVIDUAL THROW KITS

- Compact vacuum-sealed pouch with essential point of wounding equipment provided for rapid deployment to on-scene responders during a Mass Casualty Incident
- Designed to treat bleeding from penetrating and other traumatic injuries
- Small cube space fits into a BDU pocket, vest pouch, or individual bag

SPECIFICATIONS

H 7 in. x W 5 in. x D 2.5 in.
Weight: 8 oz

INDIVIDUAL THROW KITS:



INDIVIDUAL THROW KITS

OPTIONS	ITEM #
Basic	80-1104
with Combat Gauze LE	80-1105

EX-TRACT STRAP

Point of Wounding Extraction

- Point of wounding extraction
- Back pack carry
- Hasty harness for traditional and confined space lowering and hauling



SPECS

Packaged:
H 18 in.
x W 3.5 in.
x D 1 in.

Extended:
L 58 in.

Weight: 9.3 oz

EX-TRACT STRAP

ITEM #
90-0082

MCI RAPID RESPONSE THROW KITS

Hands Free, Point of Wounding Solution with the enhanced ability of treating a large number of injured patients in a large mass casualty incident

- Rapid response sling bag designed to take the guesswork out of the management of a critical incident by providing the rapid distribution of necessary life saving equipment over a large scene
- Includes 6 Individual Throw Kits with essential life saving bleeding control equipment designed for rapid deployment to on-scene responders
- Nylon sling bag features exterior MOLLE-style attachments so the user can customize the bag according to their needs
- Provides ample storage for extra gear such as a Rescue Harness for casualty evacuation (sold separately)

SPECIFICATIONS

H 12 in. x W 12 in. x D 6 in.
Shoulder strap: 28 in. to 50 in.
Weight: 3.5 lb

MCI RAPID RESPONSE THROW KITS:

OPTIONS	ITEM #
Basic	80-1106
with Combat Gauze LE	80-1107



Medical Kits

SHOULDER IFAK

- Essential life-saving tools to address injuries caused by penetrating or blast trauma
- Compact hands-free individual first aid kit providing Basic or Advanced Life Support medical gear
- Easily identifiable rugged nylon bag with shoulder strap



SPECIFICATIONS


Bag: H 9 in. x W 8.25 in. x D 2 in.
Strap: L 17 in. x W 1 in.
Weight: 1.55 lb

SHOULDER IFAK

LEVEL	ITEM #
Basic	80-1124
Combat Gauze LE	80-1125
Advanced	80-1126*


*Advanced Kits are R_X ONLY

Every Shoulder IFAK includes these core items:


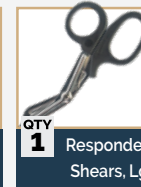


COMBAT GAUZE

also includes:

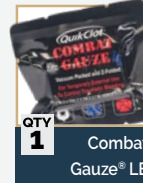

ADVANCED also includes:

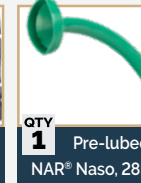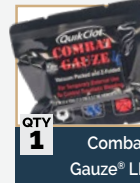

RAPID RESPONSE KIT RTF Edition

- Equipped to provide initial triage of a mass casualty incident and provide treatment for 4 or more patients
- Includes 4 C-A-T® tourniquets with holders and nylon dispenser with 4 rolls of color-coded triage tape

RAPID RESPONSE -RTF

COLOR	ITEM #
● BLK	85-0716
● RED	85-0717

SPECIFICATIONS

H 12 in. x W 12 in. x D 6 in.
Shoulder strap: 28 in. to 50 in.
Weight: 6 lb 4.5 oz

RAPID RESPONSE KIT - RTF EDITION:



Red Rapid Response RTF Kit comes standard with orange C-A-T® Tourniquets, red holders and a red tape dispenser

Triage Tape, Set of 4 (Item# ZZ-0591) also available at www.NARescue.com



RAID KIT V2



- RAID V2 pouch's unique take on the industry standard hanger combined with the utility of a fanny pack allows user to decide how to wear it based on their needs during self-aid/buddy-aid
- Designed to meet the need for an individual first aid kit with multiple carry options
- First Spear Tube makes reattaching the pouch as easy as removal



RAID KIT V2

COLOR	OPTION	ITEM#
BLK	Basic	80-1545
COY	Basic	80-1546
MTC	Basic	80-1547
RGR	Basic	80-1548
BLK	Basic w/ CG	80-1549
COY	Basic w/ CG	80-1550
MTC	Basic w/ CG	80-1551
RGR	Basic w/ CG	80-1552
BLK	Basic w/ CELOX™	80-1625
COY	Basic w/ CELOX™	80-1626
MTC	Basic w/ CELOX™	80-1627
RGR	Basic w/ CELOX™	80-1628
BLK	Advanced	80-1553*
COY	Advanced	80-1554*
MTC	Advanced	80-1555*
RGR	Advanced	80-1556*
BLK	Adv. w/ CELOX™	80-1629*
COY	Adv. w/ CELOX™	80-1630*
MTC	Adv. w/ CELOX™	80-1631*
RGR	Adv. w/ CELOX™	80-1632*

Every RAID V2 Kit includes these core items:



PLUS+

BASIC also includes:



BASIC with COMBAT GAUZE®

also includes:



BASIC with CELOX™

also includes:




ADV. with COMBAT GAUZE®

also includes:



ADV. with CELOX™

also includes:

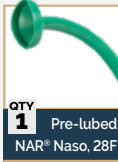


SPECIFICATIONS

H 8 in. x W 7 in. x D 2 in.
Weight: 3 lbs 3.2 oz.

MINI MEDIC KIT

- Compact, multi-purpose trauma aid bag with vital medical equipment for point of wounding treatment of multiple injuries
- Two external pockets with two external tourniquet holders
- One external pouch with additional drop-down pouch for use with breaching tools or storage for a hydration pack
- Expansion zipper to increase load capacity
- External MOLLE on sides and bottom of bag / Internal MOLLE rows on inside back of bag
- Three clear tear away pouches on the inside main compartment
- Two internal mesh pouches and one larger internal pouch on inner lid
- Internal tape roll strap



Every **MINI MEDIC KIT** includes these core items:

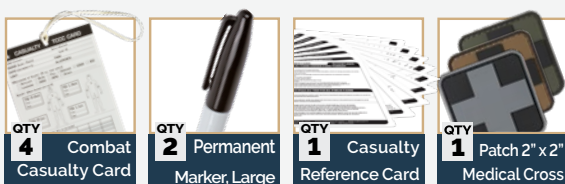



**PLUS+
BASIC**
also includes:



**COMBAT
GAUZE**
also includes:

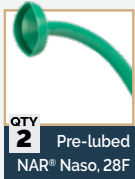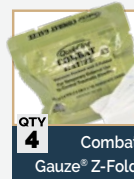

ADVANCED
also includes:

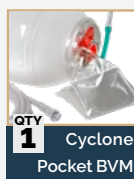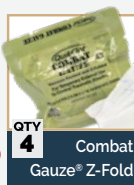

Bag only option
also available at
NARescue.com

SPECIFICATIONS

H 14 in. x
W 11 in. x D 5 in.
(Depth expandable
to 7.25 in.)
Weight: 12.5 lb
13 lb (adv.)

Inside Back
Compartment for
stowing extra
equipment

Breaching
Equipment
Sold
Separately

MINI MEDIC KIT

COLOR	LEVEL	ITEM #
BLK	Basic	80-0979
COY	Basic	80-0980
ODG	Basic	80-0981
RGR	Basic	80-1003
MTC	Basic	80-1231
BLK	Combat Gauze®	80-0985
COY	Combat Gauze®	80-0986
ODG	Combat Gauze®	80-0987
RGR	Combat Gauze®	80-1005
MTC	Combat Gauze®	80-1232
BLK	Advanced	80-0982*
COY	Advanced	80-0983*
ODG	Advanced	80-0984*
RGR	Advanced	80-1004*
MTC	Advanced	80-1233*

*Advanced Kits are Rx ONLY

Mass Casualty Incident **MCI WALK**

- All-in-one trauma treatment & evacuation platform for treatment of multiple trauma patients
- Choose from two different litter configurations: 8 non-rigid, disposable QuikLitters or 1 full-size, quad-fold, lightweight MedEVAC Litter
- Pocket system holds self- or buddy-aid point-of-wounding medical supplies
- Bag features back pack straps for easy transport from the vehicle to the point of need
- Ideal for use as a vehicle trauma kit

(80-0565/80-0458) Red MCI WALKs contain orange C-A-T® Tourniquets

(#80-0458) Red MCI WALK with MedEVAC Litter contains an orange MedEVAC Litter

MCI WALK

COLOR	LITTER	ITEM #
● BLK	8 QuikLitters	80-0564
● RED	8 QuikLitters	80-0565
● BLK	MedEVAC Litter (Black)	80-0457
● RED	MedEVAC Litter (Orange)	80-0458

SPECIFICATIONS

H 22.75 in. x W 16.5 in. x D 11.5 in.
Weight with 8 QuikLitters: 18 lb 12 oz
Weight with MedEVAC Litter: 29 lb 15 oz

MCI WALK with QUIKLITTERS:

 QTY 8 C-A-T® Tourniquet	 QTY 8 ETD® 6 in. Dressing	 QTY 5 Combat Gauze® LE	 QTY 8 NAR® Gauze, Compressed	 QTY 5 Pre-lubed NAR® Naso, 28F
 QTY 8 HyFin® Vent. Twin Pk	 QTY 8 NAR® ELITE Litter	 PAIR 8 Nitrile Gloves, Lg	 QTY 20 Combat Casualty Card	 QTY 5 Permanent Marker, Large

MCI WALK with MEDEVAC LITTER:

 QTY 8 C-A-T® Tourniquet	 QTY 8 ETD® 6 in. Dressing	 QTY 5 Combat Gauze® LE	 QTY 8 NAR® Gauze, Compressed	 QTY 5 Pre-lubed NAR® Naso, 28F
 QTY 8 HyFin® Vent. Twin Pk	 QTY 1 Rigid Litter MedEVAC	 PAIR 8 Nitrile Gloves, Lg	 QTY 20 Combat Casualty Card	 QTY 5 Permanent Marker, Large

Bag only option also
available at NARescue.com



online video
available

< scan to view

Medical Kits

MCI-TKD

Mass Casualty Incident
Throw Kit Dispenser

- Contains an inner main compartment specifically designed for quick deployment of up to eight (8) point of injury kits
- One-handed quick release, silent zippers
- G-Hook backpack straps can switch from 2 shoulder straps to bandolier style.
- Easy backpack carry for mass casualty events
- Quickly deploy bleeding control kits in emergency situations
- Side pockets have C-A-T® Tourniquets and Emergency Trauma Dressings (ETD®) for maximum preparedness for multiple injuries
- Front pockets for Litter, Breathing, and wound packing scenarios



Bag only option
also available at
NARescue.com

Every MCI-TKD includes these core items:



PLUS+

BASIC
also includes:



w/BCD
also includes:



w/CELOX™
also includes:



MCI-TKD Bleeding Control Throw Kit includes:



PLUS+

BASIC
also includes:



w/BCD
also includes:



w/CELOX™
also includes:



SPECIFICATIONS

H 21 in. x W 13 in. x D 8 in.
Weight: 15 lbs.

MCI-TKD

COLOR	LEVEL	ITEM #
● BLK	Basic	80-1749
● RED	Basic	80-1750
● BLK	w/ BCD	80-1751
● RED	w/ BCD	80-1752
● BLK	w/ CELOX™	80-1753
● RED	w/ CELOX™	80-1754



C-A-T® Combat Application Tourniquet Gen 7

- **Single routing buckle** simplifies effective slack removal and allows for fast application which decreases blood loss
- **Applied same way** to both upper and lower extremities
- True **one-handed tourniquet**
- **Scientifically proven** 100% effective in occluding blood flow in both upper and lower extremities by the USAISR and recommended by CoTCCC
- NEW welding process for **Windlass Strap** features strong adherence of strap to clip (Released 2025)



online video available
< scan to view



The **FASTEST, SAFEST, Most EFFECTIVE** Prehospital Field Tourniquet!



COMBAT APPLICATION TOURNIQUET

COLOR	STYLE	ITEM #	NSN #
● BLK	Tactical	30-0001	6515-01-521-7976
● ORG	Rescue	30-0023	6515-01-612-7129
● BLU	Trainer	30-0033	6910-01-560-2972

COMBAT GAUZE® LE

- Promotes clotting within minutes, allowing responders to control bleeding in an efficient and timely manner.
- Sterile, non-woven Z-fold gauze impregnated with kaolin for temporary external control of bleeding
- Features an embedded x-ray detectable strip
- Produces NO Heat
- 8 tear points for quicker application



SPECS

Gauze: L 4 yd x W 3 in.
Pkg: H 5 in. x W 5.5 in. x D 0.625 in.
Weight: 0.8 oz

COMBAT GAUZE® LE

ITEM

30-0109

NAR® S-ROLLED GAUZE

- Unique dispenser allows controlled linear feed application
- Used as backing gauze for hemostatic agents, to control hemorrhage
- 6-ply premium 100% cotton



SPECIFICATIONS

Gauze: W 4.5 in. x L 4.1 yd
Packaged: L 3 in. x W 2.5 in. x D 1 in.
Weight: 1.2 oz

S-ROLLED

ITEM

30-0003

NAR® COMPRESSED GAUZE

- 6-ply, premium rolled 100% cotton
- Used as backing gauze for hemostatic agents, to control hemorrhage
- Ultra-compact cube space



SPECIFICATIONS

Gauze: W 4.5 in. x L 4.1 yd
Packaged: L 3 in. x W 2 in. x D 1 in.
Weight: 1.2 oz

COMPRESSED

STYLE	ITEM #
Responder	30-0182
Standard	30-0052

CELOX™ RAPID Hemostatic Gauze



- Stops severe arterial bleeding with only one minute of compression applied
- Compact Z-folded format making it easy to use & carry
- Generates no heat

SPECIFICATIONS

Packaged:
H 5.25 in. x W 4.25 in. x D 0.5 in.
Gauze: W 3 in. x 5 ft
Weight: 1.1 oz

CELOX™ RAPID

ITEM

30-0069

NAR® WOUND PACKING GAUZE

- Ultra compressed, Z-fold gauze indicated for wound packing
- Soft, white 1-ply gauze with an embedded x-ray detectable line
- Z-fold configuration for controlled application



SPECIFICATIONS

Gauze: W 4.5 in. x L 4.1 yd
Packaged: L 3 in. x W 2 in. x D 1 in.
Weight: 1.2 oz

WOUND PACKING

ITEM

30-0054

NAR® RESPONDER ETD® Emergency Trauma Dressing

- Sterile pressure bandages that deliver extremely effective pressure and excellent absorption
- Suitable for applying immediate direct pressure
- Easy-to-use hook & loop securing device (30-0179 also features C-Clip)
- Vacuum-sealed, low-cube packaging

RESPONDER EMERGENCY TRAUMA DRESSINGS

STYLE	BDG.	PAD	WGT	PACKAGE	ITEM #
4in.Flat	4 x 60	4x6	2 oz	L 4.25 in. x W 2.75 in. x D 0.88 in.	30-0171
4in.(Mini)	4 x 42	4x6	2.4 oz	L 5 in. x W 3.5 in. x D 0.75 in.	30-0179
4in.Roll	4 x 60	4x6	2 oz	L 4.5 in. x W 2.5 in. x D 1.75 in.	30-0167
6in. Flat	6 x 56	6x6	3.31 oz	L 6.25 in. x W 3.75 in. x D 0.88 in.	30-0172
6in. Roll	6 x 60	6x6	2.9 oz	L 6 in. x W 2.5 in. x D 2.1 in.	30-0168

NAR® PRE-LUBRICATED NASOPHARYNGEAL AIRWAY 28F

- **NEW ITEM:** Pre-lubricated* Nasopharyngeal Airway offers for more efficient insertion and decreased procedure prep time
* Water based lubricant (-5°F freezing point)
- NAR® Pre-lubricated Nasopharyngeal airways are **MADE IN USA!**
- Nasopharyngeal airways are designed to be more compatible with human anatomy than the standard oropharyngeal airway
- Rounded bevel tip allows for gentle, optimal insertion
- Soft material avoids damage to nasal passage
- Latex-free and sterile

HYFIN® VENT COMPACT Twin Pack

- Two Chest Seals packaged approximately 25% smaller than our standard HyFin®
- 3-channel pressure relief vents
- Two Chest Seals designed for the treatment of entry/exit or multiple penetrating wounds
- Large red tab designed for single step, peel and apply application and for burping the wound if necessary

HYFIN® VENT COMPACT

ITEM

10-0042

SPECIFICATIONS

Packaged:
H 3.88 in x W 6.6 in x D 0.13 in.
Device (excludes tab):
4.75 in. x 4.75 in.
Weight: 1.55 oz



SPECIFICATIONS

L 6 in. x W 16 in. x D 0.5 in.
Weight: 0.5 oz

NAR® Pre-Lubricated Nasopharyngeal Airways

OPTION	ITEM#	NSN#
Single	10-0060	6515-01-710-9764
Box of 10	10-0067	6515-01-710-9783



Other sizes (20F-36F) and styles available at NARescue.com

NAR® SURVIVAL BLANKET

- Lightweight, durable space blanket made of a combination of mylar and polyethylene, perfect for regulating body temperature
- Reflective silver color helps to put radiant heat back into the body
- Wind and waterproof, further aiding in retaining warmth
- Full body coverage



NAR® SURVIVAL BLANKET

ITEM #	NSN #
ZZ-0038	4240-01-697-8869

SPECIFICATIONS

Packaged: L 5 in. x W 3 in.
Blanket: L 84 in. x W 56 in.
Weight: 1.84 oz.

Responder NITRILE GLOVES

- Form-fitting, ambidextrous design with textured fingertips
- Thick palm ensures durability
- Textured for improved wet and dry gripping
- Longer and thicker than standard nitrile gloves
- Also available in box of 100 at www.NARescue.com

RESPONDER NITRILE GLOVES

SIZE	PKG QTY	ITEM #
S	Bag 25 pair	ZZ-1059
M	Bag 25 pair	ZZ-1060
L	Bag 25 pair	ZZ-1043
XL	Bag 25 pair	ZZ-1061

SPECIFICATIONS

Bag (25 pair): L 4.75 in. x W 6 in. x D 4 in.
Weight: 1 lb 1 oz



Patient Evacuation | Litters & Stretchers

MEDEVAC LITTER

- Weighs 13.5 lb and can carry a payload up to 1,200 lb
- Built-in patient securing buckles
- Has six built-in IV attachment points
- Ergonomically designed handles
- Patented joint configuration locks the Leg Crossbar in place and prevents unintentional collapsing
- 400D/PU nylon litter bed



SPECIFICATIONS

Folded: H 21.5 in. x W 8.5 in. x D 7.5 in.
Open: H 84.5 in. x W 22.5 in. x D 4 in.
Weight: 13.5 lbs.
Static Payload: 1,200 lbs.
Ground Clearance: 3.9 in.

MEDEVAC

COLOR	ITEM #
BLK	60-0051
ORG	60-0052



NAR® ELITE RESPONDER LITTERS

Evacuation Litter for Incident Tactical Extraction

- Multi-purpose, low cost disposable evacuation platforms for transport & transfer
- Water-resistant lightweight polyethylene construction with 10 large hand-holds
- Super compact LITE is made of nylon & vacuum packed to fit in a cargo pocket
- The XL features 12 hand-holds and 2 built-in retention straps
- ELITE - RED features signature red nylon webbing & stitching



DISPOSABLE QUIKLITTERS

MODEL	LITTER BED	WGT	PAYLOAD	ITEM #
LITE*	L 74 in. x W 27 in.	13 oz	750 lb	60-0071
Standard	L 76 in. x W 30 in.	17 oz	500 lb	60-0061
RED	L 76 in. x W 30 in.	17 oz	500 lb	60-0080
XL	L 80 in. x W 40 in.	2.2 lb	1000 lb	60-0065

*Nylon Construction

First responders have increasingly become **the target of violence** in recent times.

Nationwide, **violent crime is rising** at an alarming rate, translating to a higher call volume. With the **sharp decline** in the **police presence** which has historically provided scene safety for Fire/EMS, First Responders are experiencing a **dramatic increase** in **dangerous** and **potentially violent situations** in their everyday operations. Unfortunately rescuers in the Fire/EMS arena can be perceived as just another target for potential life threatening violence along with their law enforcement counterparts.

Now more than ever, your **personal safety is paramount** in the performance of your duties. Wearing **ballistic protective equipment (BPE)** is one of the best ways **to help protect yourself** in a violent environment.

PROTECT Yourself & SAVE Lives...

Law enforcement transitioned to wearing body armor on a daily basis in the 1980s because they realized it was an essential component of their daily equipment.

It saves lives! In the same way, BPE is critical for Rescue Task Force (RTF) and Tactical Emergency Medical Service (TEMS). Fire/EMS should now consider Ballistic Protection Equipment for daily wear on any call that could result in violence or an unsafe environment.

Protecting Yourself means you can save others AND live to serve another day!

Sign up to get updates and special offers when you...

#JOIN THE MISSION



NOTICE

Ballistic protective equipment (BPE) products are controlled by the Export Administration Regulations (EAR). It may not be exported, sold, transferred, transported, or shipped outside of the U.S., either in its original form or after being incorporated into other items, without first obtaining approval from the U.S. government or as authorized by U.S. law and regulations.



Explore all your options at

NORTH AMERICAN RESCUE®

www.NARescue.com • 888.689.6277