

MASSIVE HEMORRHAGE | AIRWAY | RESPIRATION | CIRCULATION | HYPOTHERMIA

EMS *Catalog*





NORTH AMERICAN RESCUE®

www.NARescue.com • 888.689.6277

Since 1996, North American Rescue® has focused on ***increasing survivability*** on both the battlefield and the homefront.

With our innovative medical/rescue equipment, NAR® continues to lead the way by delivering the ***tailored, mission-critical healthcare solutions*** needed by today's first responders to counter the rapidly evolving threats of the 21st century. See what a lifesaving difference our combat proven ***Products with a Mission®*** make in the prehospital environment!

PRODUCTS WITH A MISSION®



BE THE DIFFERENCE®

www.NARescue.com





We are dedicated to decreasing preventable death by providing our products at special pricing to our EMS and pre-hospital life savers. To see if you qualify for this discount, please contact us or fill out our online NAR® Pricing Application. Our customer service department is standing by to help!

888.689.6277



M.A.R.C.H.

Massive Hemorrhage **4-6**

Airway Management **7**

Respiratory Support **8-11**

Circulation & Vascular Access **11-15**

Hypothermia Management **16-17**

Trauma Kits

Custom Kitting Solutions **18-19**

IFAKS **20-21**

Multiple Injury Kits **22-24**

Extras

Splinting & Evacuation **25**

Treatment Support Accessories **26-27**

Education & Training **28-29**

Massive Hemorrhage | Tourniquets

C-A-T® Combat Application Tourniquet Gen 7

- **Single routing buckle** simplifies effective slack removal and allows for fast application which decreases blood loss
- **Applied same way** to both upper and lower extremities
- True **one-handed tourniquet**
- **Scientifically proven** 100% effective in occluding blood flow in both upper and lower extremities by the USAISR and recommended by CoTCCC
- NEW welding process for **Windlass Strap** features strong adherence of strap to clip (Released 2025)



online video available

< scan to view

The FASTEST, SAFEST, Most EFFECTIVE Prehospital Field Tourniquet!

COMBAT APPLICATION TOURNIQUET

COLOR	STYLE	ITEM #	NSN #
Black	Tactical	30-0001	6515-01-521-7976
Orange	Rescue	30-0023	6515-01-612-7129
Blue	Trainer	30-0033	6910-01-560-2972



SPECIFICATIONS

Packaged: H 6.5 in. x W 2.4 in. x D 1.5 in.
Device: L 37.5 in.
Weight: 2.7 oz.

RIGID C-A-T® TQ CASE *with Shirt Shield*

- Allows carriage of a Gen 7 C-A-T® Tourniquet on your duty/patrol/pistol belt (Tourniquets sold separately)
- Shirt Shield protects uniform shirts from undue friction wear
- Blade-Tech TekLok Belt Attachment fits up to a 2.25 in. belt
- Can be mounted either vertically or horizontally, or to a Safariland Drop-leg Shroud
- Injection Molded Nylon Case stands up to rugged military usage
- Cross medical identifier patch included

RIGID C-A-T TQ CASE

COLOR	ITEM #
Black	80-0944
Black, Basket Weave	80-0962
Black, Polished	80-0963
Coyote	80-0945
Ranger Green	80-0946
Navy	80-1161

SPECIFICATIONS

H 5.88 in. x W 2.5 in. x D 2 in.
Weight: 5.36 oz

Explore other TQ carrier options available at **NARescue.com**

RIGID C-A-T® ITW TQ CASE

Inside-the-Waistband Tourniquet Case

- Holds the Gen 7 C-A-T® Tourniquet (Sold Separately)
- Designed with built-in Shirt Shield
- Drainage hole
- 1.5" and 1.75" IWB Overhook Mount
- Designed, manufactured, and assembled in the USA.

RIGID CAT® TQ CASE

COLOR	ITEM #
BLK	80-1488
GRY	80-1489

SPECIFICATIONS

Packaged: L 5 in. x W 3.25 in. x D 1.625 in.
Weight: 3.0 oz

Rigid Case MOLLE ATTACHMENT

- 1 Malice Clip
- 1 MOLLE Weaveback
- 2 Barrel Nuts
- 2 Screws

RIGID MOLLE ATTACHMENT

ITEM #
ZZ-1696

SPECIFICATIONS

H 4 in. x W 1 in. x D 0.5 in.
Weight: 0.4 oz.



LEATHER C-A-T® HOLDERS



- Designed to meet the specific duty gear needs of law enforcement officers - for use on daily duty belt to fit within organizational uniform standards
- Securely fastens to your duty belt with easy one-handed access to allow you to take control of extreme medical situations in seconds
- Constructed of black full grain leather in different finishes, attachment options, and snaps colors
- High quality snaps ensure secure placement on your duty belt and easy accessibility to your tourniquet
- Ideal for the C-A-T® Gen 7 but also accommodates other popular windlass TQs
- Attachment options: Clip, Loop or Hook & Loop Strap attachment on select finishes and snap colors



UNIFORM DUTY BELT LEATHER TQ HOLDER

FINISH	ATTACHMENT	SNAP	ITEM #
High Gloss	Clip	Silver	30-0194
Basket Weave	Loop	Silver	30-0195
Full Grain	Loop	Silver	30-0197
Basket Weave	Strap Attach	Silver	30-0228
High Gloss	Strap Attach	Silver	30-0229

SPECIFICATIONS

L 5.75 in. x W 2 in. x D 2 in.
Weight: 3.2 oz.

JETT® Junctional Emergency Treatment Tool

- Immediately treat high extremity wounds in the thigh/groin area that is not amenable to TQ application
- Simultaneous occlusion of blood flow to both lower limbs that will not impede respiration
- Patent 9,492,177



SPECIFICATIONS

R_X ONLY

Packaged:
H 5.5 in. x W 7 in. x D 3.75 in.
Overall Device: L 56.75 in.
Buttocks Pad: L 19 in. x W 6.5 in.
Compression Pad:
L 4.85 in. x W 2.5 in. x D 11 in.
T-Handle:
H 4.62 in. x W 4.04 in. x D 0.63 in.
Weight: 1 lb 9.6 oz

JETT

ITEM #	NSN #
30-0088	6515-01-616-5841

online video
available
scan to view >



Massive Hemorrhage | Bandages & Gauze

NAR®

COMPRESSED GAUZE

- Low-cube, ultra-compact sterile gauze
- Premium rolled 100% cotton gauze for maximum absorption



COMPRESSED GAUZE

STERILE RED TIP

STYLE	ITEM #	NSN #
Responder	30-0182	----
Standard	30-0052	6510-01-503-2117

SPECS

Gauze: W 4.5 in. x L 4.1 yd
Packaged: L 3 in. x W 2.5 in. x D 1 in.
Weight: 12 oz

NAR®

WOUND PACKING GAUZE

- Ultra compressed, Z-fold gauze indicated for wound packing
- Soft, white 1-ply gauze with an embedded x-ray detectable line



SPECS

STERILE RED TIP

Gauze: W 3 in. x L 5 yd
Packaged: H 3 in. x W 2 in. x D 1 in.
Weight: 12 oz

WOUND PACKING GAUZE

ITEM #	NSN #
30-0054	6510-01-695-5824

NAR®

RESPONDER ETD®

Emergency Trauma Dressing

- Sterile pressure bandages that deliver extremely effective pressure & excellent absorption
- Suitable for applying immediate direct pressure
- Equipped with a non-adherent pad
- Rugged vacuum-sealed, low-cube packaging with easy-open tear notches
- Quick grip roll control system prevents accidental unrolling
- Easy-to-use hook & loop securing device (30-0179 also features C-Clip)



RESPONDER ETD®S

STERILE RED TIP

STYLE	BDG.	PAD	WGT	PACKAGE	ITEM #
4in.Flat	4 x 60	4x6	2 oz	L 4.25 in. x W 2.75 in. x D 0.88 in.	30-0171
4in.(Mini)	4 x 42	4x6	2.4 oz	L 5 in. x W 3.5 in. x D 0.75 in.	30-0179
4in.Roll	4 x 60	4x6	2 oz	L 4.5 in. x W 2.5 in. x D 1.75 in.	30-0167
6in. Flat	6 x 56	6x6	3.31 oz	L 6.25 in. x W 3.75 in. x D 0.88 in.	30-0172
6in. Roll	6 x 60	6x6	2.9 oz	L 6 in. x W 2.5 in. x D 2.1 in.	30-0168

Massive Hemorrhage | Hemostatics

QUIKCLOT® EMS | COMBAT GAUZE® | COMBAT GAUZE® LE

- Sterile, non-woven, gauze impregnated with kaolin for temporary external control of bleeding
- Embedded x-ray detectable strip
- Combat Gauze® products are Z-Folded



SPECIFICATIONS

STERILE

Standard & LE Gauze:

L 4 yd x W 3 in.
Pkg: H 5 in. x W 5.5 in. x D 0.625 in.
Weight: 0.8 oz

QuikClot EMS *Rx ONLY

L 4 in. x W 4 in.
Pkg: H 5.5 in. x W 5 in. x D 0.2 in.
Weight: 0.8 oz

QuikClot EMS (Rolled):

L 4 ft. x W 3 in.



HEMOSTATICS

VERSION	ITEM #	NSN #
Combat Gauze LE	30-0109	----
Combat Gauze	30-0039	6010-01-562-3325
Combat Gauze XL	30-0070	6510-01-600-4705
QuikClot EMS*	30-0095	----
QuikClot EMS (Rolled)	30-0166	

BAC-PACK™ Bougie-Aided Cricothyrotomy Pack

- All-in-one cricothyroidotomy kit for use in any environment
- Pre-inserted bougie in tube allows for a more rapid, easy procedure
- Tube meets TCCC recommendations for bougie-aided open surgical technique
- #10 Scalpel included for safety and ease of use
- Small, durable pouch keeps the items organized and ready to use
- Lightweight and ruggedized
- Components configured by order of use
- Scalpel, tube, syringe all provided sterile
- Patent pending

BAC-PACK



SPECIFICATIONS


Packaged: L 5.5 in. x H 3.625 in. x D 0.5 in. Tube length: 67.0mm
 Inner tube diameter: 6.0mm Weight: 0.84 oz
 Outer tube diameter: 8.0mm US Patents: Pending

BAC-PACK

ITEM

10-0052

NAR® Cricothyrotomy Kit with BOUGIE INTRODUCER

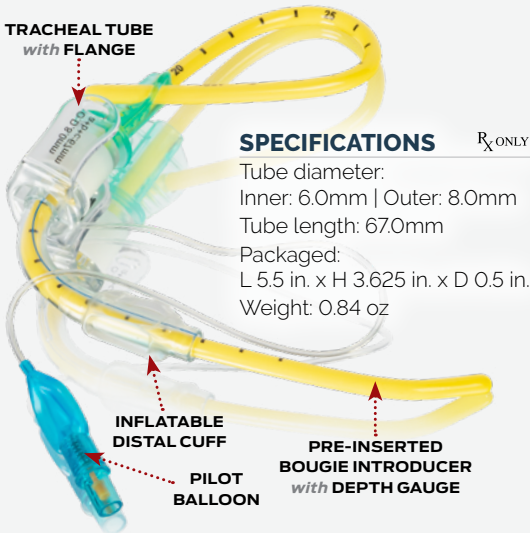

- Convenient device that gives the user a pre-inserted bougie for introduction of the tube into the trachea and indication of correct placement prior to insertion
- Adjustable for maximum user preference
- Latex free tube features 6.0mm inner diameter and inflatable cuff
- Contains the same preferred flexible tube as NAR®'s Tactical Crickit®

CRICOTHYROTOMY KIT


ITEM

10-0048

NAR® NASOPHARYNGEAL AIRWAY 28F (117mm)



- **NEW ITEM:** Pre-lubricated for more efficient insertion
- Nasopharyngeal airways are designed to be more compatible with human anatomy than the standard oropharyngeal airway
- Other sizes (20F-36F) and styles available at NARescue.com

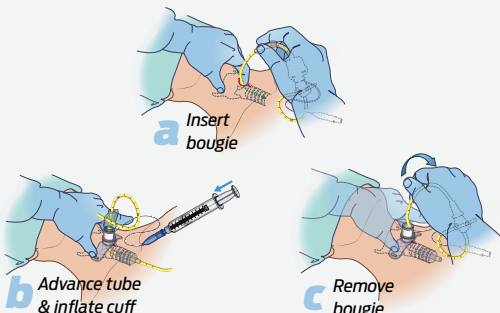

NAR® Nasopharyngeal Airway, 28F (117mm)

OPTION	ITEM #	NSN #
Pre-lubricated	10-0060	6515-01-710-9764
Non-lubricated	10-0061	6515-01-710-9815
Pre-lubricated Box	10-0067	6515-01-710-9783
Non-lubricated Box	10-0068	6515-01-710-9848

SPECIFICATIONS

Packaged Device:
 L 7.5 in. x H 3.5 in.
 Device Weight: 0.02 lb.
 Box: L 9 in. x W 4 in. x D 4 in.
 Box Weight: 0.5 lb.




HYFIN® CHEST SEALS



HYFIN® VENT COMPACT, *Twin Pack*

- Two Chest Seals packaged approximately 25% smaller than our standard HyFin®
- 3-channel pressure relief vents
- Two Chest Seals designed for the treatment of entry/exit or multiple penetrating wounds
- Large red tab designed for single step, peel and apply application and for burping the wound if necessary

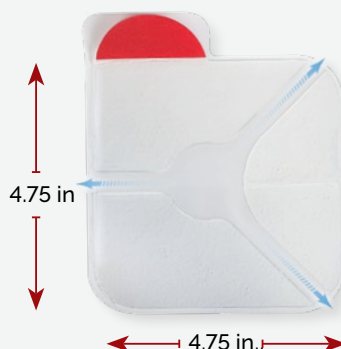

HYFIN VENT COMPACT

ITEM #
10-0042

SPECIFICATIONS

Device (excludes tab): 4.75 in. x 4.75 in.
Weight: 155 oz
Packaged: H 3.88 in. x W 6.6 in. x D 0.13 in.



online video available

scan to view >



HYFIN® VENT *Twin Pack*

- Two individually packaged, vented/valve-type chest seals in one unique package designed for treatment of both entry/exit wounds or multiple penetrating injuries
- Patented 3 channel pressure relief vents
- Ruggedized dual perforated pouch designed to allow the provider to open only one dressing at a time, as needed
- Advanced adhesive technology designed to provide a superior seal

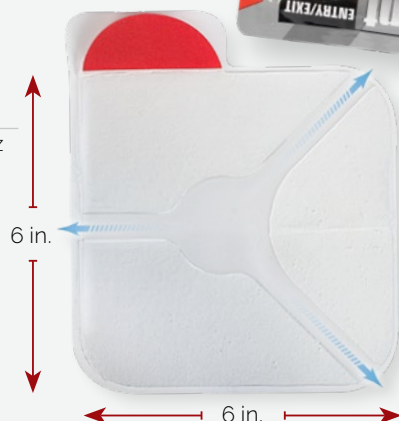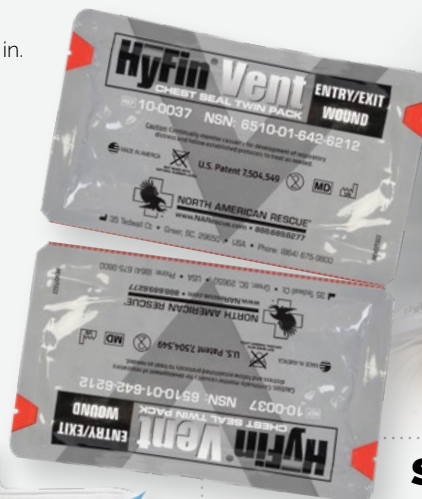

SPECIFICATIONS

Device (excludes tab): 6 in. x 6 in. Wgt: 2.5 oz
Packaged: H 7.5 in. x W 4.5 in. x D 0.25 in.
Unfolded: H 7.5 in. x W 9 in. x D 0.13 in.

HYFIN VENT TWIN PACK

ITEM #	NSN #
10-0037	6510-01-642-6212

SOTAIR CLINICAL DEVICE 6 Pk

- Fits between a manual resuscitator and a patient mask/airway
- Regulation of airflow minimizes the risk of over-pressurization, over-ventilation and hyperventilation



SPECIFICATIONS

Diameter: L 6.18 in. x W 4.29 in. x D 2 in.
Weight: 0.35 lbs.

SOTAIR - 6PK

ITEM #
10-0083



Includes 6 ft 6 in. of oxygen tubing & rugged, hard plastic, twist-top case

CYCLONE® Improved POCKET BVM

- Unique collapsible design packaged 25% smaller than conventional BVMs
- Quick to deploy and does not need to be manually inflated prior to use

SPECIFICATIONS

Tidal (Balloon) Volume: 1600 ml
Diameter: 5.3 in. x H 2.9 in.
Weight: 17.4 oz

CYCLONE® POCKET BVM

ITEM #	NSN #
10-0019	6515-01-568-0193

ENHANCED ARS™ Air Release System

3.25 inch



10 ga ENHANCED ARS™

Needle Hub
(no flash chamber)

10ga fenestrated catheters were shown to offer marked increase in air outflow and a corresponding decrease in time to relief of tension, when compared to all NDC Catheters

Centimeter Markings
(with contrasting colors that are NVG compatible) for safer placement utilizing "thoracic depth control"

Modified Needle Tip
"scalpel" sharp, bi-bevel, tapered needle tip, which improves penetration control through tissue of varied type and thickness

Catheter Hub
color identifies needle catheter gauge

10ga (or 14ga) 3.25-inch Flexible Catheter to decrease incidence of catheter kinking, in concert with modified insertion recommendations from leading experts in pulmonology, critical care, emergency medicine, and trauma surgery.



Fenestrations showing statistically significant improvement in tension relief and the prevention of distal catheter occlusion when compared to traditional non-fenestrated catheters



"It's time to rethink needle decompression techniques and device characteristics to optimize success"

SPECIFICATIONS

STERILE R_X ONLY

Pkg: H 5.7 in. x W 0.7 in. x D 0.8 in.
Weight: 0.6 oz
US Patents: D584410-S

ENHANCED ARS

GAUGE	ITEM #	NSN #
10 ga	10-0064	----
14 ga	10-0063	----




↑ scan to learn more online

Respiratory Support | Tension Pneumothorax

10 ga x 3.75 inch **SPEAR®** Simplified Pneumothorax Emergency Air Release

- Latest innovation in Needle Decompression technology
- Designed based on known gaps in the management of tension pneumothorax
- Enhanced needle tip for improved insertion control
- User friendly features like spin-lock connection, depth markings, and detachable low-pressure one-way valve
- Fenestrations along catheter for additional air flow
- Familiar, rugged, pen-style case for ease of storage and deployment
- Allows for either lateral decompression or anterior needle thoracostomy
- Catheter made of flexible kink-resistant material with a soft tip
- Single use only. Supplied sterile

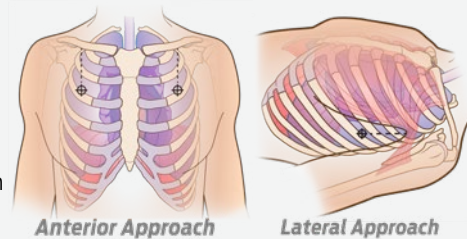

SPEAR

GAUGE	ITEM #	NSN #
10 ga	10-0051	6515-01-699-4360


Pneumothorax Decompression Indicator **CAPNOSPOT™**

- Objective, visual indication of decompression
- Real-time confirmatory color change (from blue to yellow).
- Compatible with existing decompression devices, catheters, and needles via standard Luer Lock.

SPECIFICATIONS Rx ONLY

Weight: 3.2 oz.

ITEM #	NSN #
10-0070	----

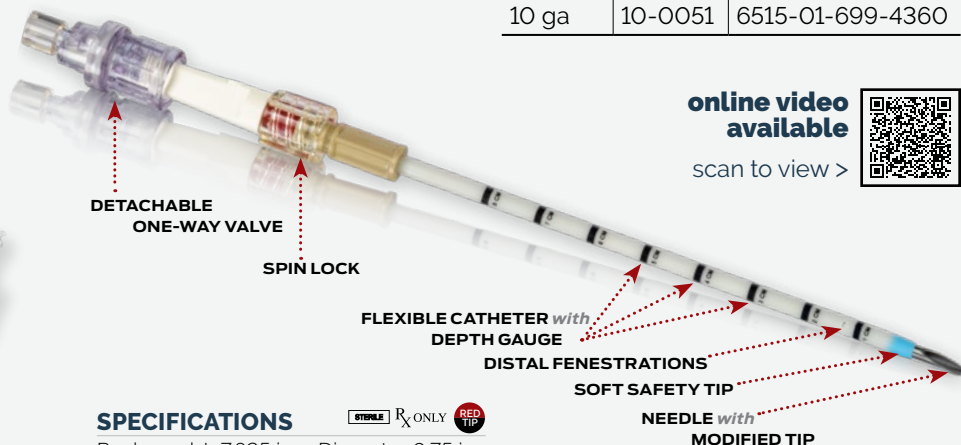

SPECIFICATIONS

Packaged: L 7.625 in. x Diameter 0.75 in.
Catheter: 10Ga x 3.75 in.
Weight: 1 oz
US Patents: 10926063, D584410-S

online video available

scan to view >



3.25 inch **ARS®** NEEDLE DECOMPRESSION KIT

- Penetrate the pleural space in 99% of casualties presenting the signs and symptoms of a tension pneumothorax
- Quality, reliable devices with the appropriate length needle/catheter
- Capless flash chamber provides immediate confirmation of correct placement
- Patented design packaged in an easy-to-open, crush-proof protective tube featuring easy-to-identify textured twist top with clip

ARS NEEDLE DECOMPRESSION KIT

GAUGE	ITEM #	NSN #
14 ga	ZZ-0056	6515-01-541-0635
10 ga	ZZ-0298	6515-01-673-1701

U.S. Design Patents: D584,409 S; D595,847S; D584,410S; U.S. Patent 7,874,426B2; © North American Rescue,® LLC., 127820, 127822, 127823, 001002372-0001, 001013940-0001, 001013940-0002



SPECIFICATIONS

Pkg: H 5.7 in. x W 0.7 in. x D 0.8 in.
Weight: 0.6 oz
US Patents: D584410-S



online video available

< scan to view



Rx ONLY

SPECIFICATIONS


Packaged: L 4.5 in. x W 7.5 in. x D 0.75 in.
Weight: 5.28 oz

NAR® SIMPLE THORACOSTOMY KIT

- Compact kit that contains all necessary equipment to perform simple thoracostomy
- Includes 3-pack skin prep for performing the procedure bilaterally if needed
- HyFin® Compact Vented Twin Pack Chest Seal to cover site if required

SIMPLE THORACOSTOMY KIT

ITEM #	UNIT
83-0027	Single

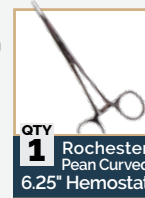

QTY 1 Rochester Pean Curved 6.25" Hemostat



QTY 1 Protected #10 Scalpel



QTY 2 4" x 4" Gauze Pad, 2pk



QTY 1 HyFin® Vent Compact, Twin




QTY 1 Chloraprep Swabsticks, 3pk



QTY 1 Permanent Marker, Small

CHEST TUBE INSERTION KIT

- Designed to provide the necessary equipment needed to perform emergency chest tube thoracostomy
- Compact and lightweight kit packaged in an easy to open, vacuum sealed pouch
- Can be stored in your kit or crash cart

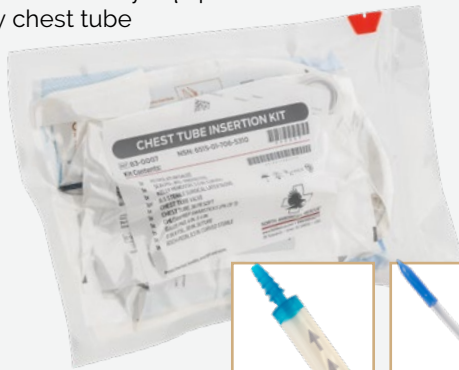

Rx ONLY

CHEST TUBE INSERTION KIT

ITEM #	NSN #
83-0007	6515-01-706-5310



QTY 1 Gauze, Petrolatum



QTY 1 Protected #10 Scalpel



QTY 1 5.5" Kelly Hemostat, Cvd



PAIR 1 Surgical Gloves, Latex



QTY 1 Valve Chest Tube



QTY 1 36 French Chest Tube




QTY 1 Chloraprep Swabsticks, 3pk



QTY 1 4" x 4" Gauze Pad, 2pk



QTY 1 O-Silk Sutures, 30 in.

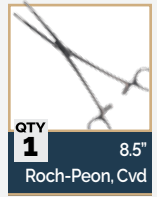

QTY 1 8.5" Roch-Peon, Cvd

SPECIFICATIONS

Packaged: L 10 in. x W 6 in. x D 2.25 in.
Weight: 11.2 oz

Circulation & Vascular Access

BOA® IV Constricting Band

- Direct peripheral IV access utilizing a revolutionary rolling action for transfusing venous blood to the distal veins while still allowing arterial flow
- Quick Click connector & two-finger Quick Release improves "time on task"
- XL version stretches approximately 50% longer than the standard BOA®



online video available

< scan to view



SPECIFICATIONS


Standard	XL	Packaged
Dimensions - L 7.5 in. to 18 in.	Dimensions - L 10 in. to 27 in.	L 4 in. x W 3 in. x D 0.9 in.
Weight - 12 oz.	Weight - 15 oz.	US Patents: 5,304,202 & 5,607,448



BOA IV CONSTRICTING BAND

SIZE	ITEM #	NSN #
Standard	30-0009	6515-01-537-2611
XL	30-0071	----

SAM[®] IO Intraosseous Access System

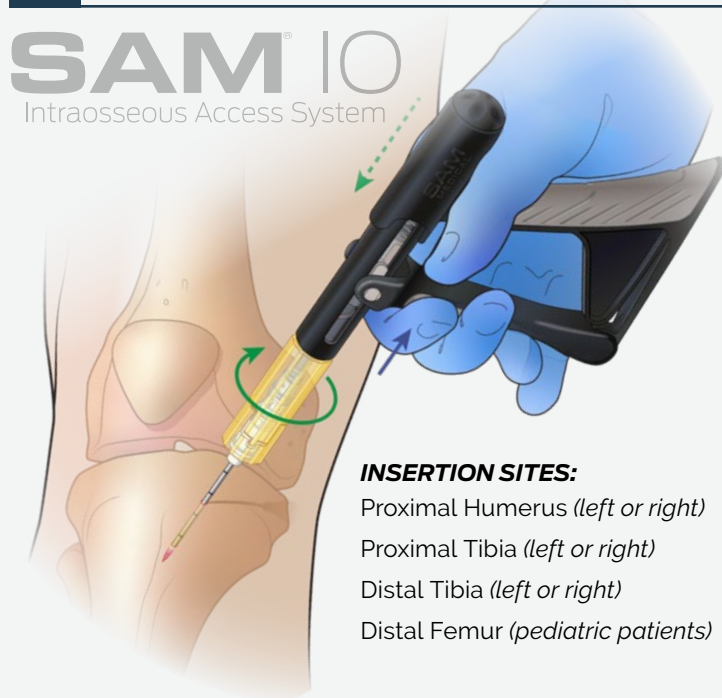

INSERTION SITES:

Proximal Humerus (left or right)

Proximal Tibia (left or right)

Distal Tibia (left or right)

Distal Femur (pediatric patients)

SAM[®] IO KIT *Intraosseous Access*

- Designed for intuitive, effective IO placement
- Includes multi-use manually operated intraosseous access driver with passive ratcheting mechanism that rotates clockwise when resistance is encountered
- Provides controlled insertion of the IO needle with a tactile feel of cortical penetration
- Includes 15 mm, 25 mm, and 45 mm Needle sets containing extension tubes and single-use sharps containment to address a broad range of patient demographics and multiple application sites
- Featuring a stabilizer with cylinder design intended to protect the catheter from dislodgement

SPECIFICATIONS R_x ONLY

H 7 in. x W 11.25 in. x D 5.5 in.

Weight: 2.32 lb

SAM IO KIT

ITEM

30-0307

SAM[®]
MEDICAL

See item# below for additional **SAM IO System** components also available for purchase:

SAM IO INTRAOSSEOUS ACCESS SYSTEM / INDIVIDUAL ITEMS

PRODUCT	UNIT	ITEM #
Driver	Single	30-0255
Needle 15mm	Single	30-0257
Needle 15mm	5-Pack	30-0258
Needle 25mm	Single	30-0259
Needle 25mm	5-Pack	30-0260
Needle 45mm	Single	30-0261
Needle 45mm	5-Pack	30-0262
Stabilizer	Single	30-0263
Stabilizer	5-Pack	30-0264
Adaptor*	Single	30-0253
Training Kit	Kit	30-0256
Storage Case	Single	30-0254

*SAM IO Adaptor is designed to accommodate alternative IO needle designs with the SAM IO Driver

SAM IO KIT



Includes **Custom Field Storage Case** with elastic loops for easy transport



Extension Tube & NeedleVISE[®]
included in each needle set

Crêdo PROMED™ SERIES FOUR CONTAINER

- Temperature-controlled portable medical transport bags that thermally protect the integrity of valuable medical supply payloads for 48 - 96 hours
- Outer Container: Ballistic nylon fabric with durable web strapping and secure buckles
- Featuring the patented Thermal Insulation Chamber (TIC™) panels with phase change material and VIP components contained within to consistently protect medical materials such as blood and platelets within 2-8°C.
- Available in 3 capacities: 2L, 4L & 8L

CRÊDO PROMED™ SERIES FOUR CONTAINER

CAPACITY	ITEM #
2L	35-0011
4L	35-0012
8L	35-0013
2L Chamber	35-0021
4L Chamber	35-0022
8L Chamber	35-0023

SPECIFICATIONS

2 Units: Volume: 2L
Duration: 48 hours
L 8 in. x W 9 in. x D 10 in.
Weight: 764 lb

4 Units: Volume: 4L
Duration: 96 hours
L 10 in. x W 10.5 in. x D 10.5 in.
Weight: 11.45 lb

8 Units: Volume: 8L
Duration: 72 hours
L 11 in. x W 15.5 in. x D 11 in.
Weight: 19.6 lb



*Ideal for use with the
Quantum™ and LifeFlow
devices to provide a complete
blood delivery package*

PAX COOLER BAG For the Crêdo ProMed™ Series Four - 2L Container

- Features the non-porous PAX Tech material that provides an easy-to-clean and decontaminate exterior
- Has a detachable section that is ideal for storing the Quantum™ Blood & Fluid Warming System and other essential IV access equipment
- Compatible with the Thermal Insulation Chamber - TIC™ and VIP components from the 2 Liter Crêdo cooler panels

PAX COOLER BAG

CAPACITY	ITEM #
2L	35-0032
4L	35-0033

SPECIFICATIONS

H 7.5 in. x W 11.5 in. x D 10 in.
Weight: 2 lbs 1.6 oz



*Crêdo ProMed™
sold separately*



Circulation & Vascular Access

APRU - 6L AUTONOMOUS PORTABLE REFRIGERATION UNIT - 6 LITER

Features

- FDA Listed & AABB Certified
- Complies with Blood Bank requirements
- Rugged design, easily transportable
- Waterproof and dustproof
- **Authorized** for use on board UH-60L, UH/HH-60M, and UH-60V Army helicopters in battery powered mode

Storage

- Internal temperature: 4.5°C ± 0.8°C
- Duration: 95 hrs @ 24°C (@ 75°F) on single battery charge
- Volume: 4 L (4.4 qt)
- Unit capacity: Six 500 ml bags

Temperature Recording

- Frequency: 4 hours
- Duration: 72 days
- Accuracy: ± 1°C, class 2 (k=2)
- Calibration: NIST Traceable

Contents not included



APRU - 6L

ITEM #

35-0034

SPECIFICATIONS

Dimensions: H 10.1 in. x W 19.3 in. x D 10.2 in.
Weight (unloaded): 29 lbs.



Contents Sold Separately

APRU TRANSPORT BOX

- Heavy duty, with laser-cut custom foam inserts
- Compatible with the Autonomous Portable Refrigeration Unit (APRU) and its accessories
- APRU and its accessories are sold separately

SPECIFICATIONS

Box Dimensions - H 30 in. x W 18 in. x D 23 in.
Weight: 41 lbs.

TRANSPORT BOX

ITEM #

35-0035

APRU VEHICLE MOUNT

- Compatible with the Autonomous Portable Refrigeration Unit (APRU) - *sold separately*
- Does not include strap



MOUNT

ITEM #

35-0036

SPECIFICATIONS

Vehicle Mount Dimensions: H 12.5 in. x W 10 in. x L 19.25 in.
Hardwire Cable Length: 15 ft.
Total Weight: 5.2 lbs.

APRU BATTERY BB2590

- Rechargeable Lithium-Ion Battery
- Compatible with the Autonomous Portable Refrigeration Unit (APRU) - *sold separately*
- Designed and manufactured in the USA
- **Nominal Voltage:** 14.4V, 28.8V
- **Ah:** In 24V Mode: 7.5 Ah, In 12V Mode: 15 Ah
- **Maximum Voltage:** 16.8V



BATT.

ITEM #

35-0039

SPECIFICATIONS

Dimensions: H 5 in. x W 2.4 in. x L 4.4 in.
Weight: 3.1 lbs.

APRU BACKPACK CARRIER WITH POUCHES

- Adds carrying capability by hand or as a backpack
- Compatible with the Autonomous Portable Refrigeration Unit (APRU) - *sold separately*
- Includes two extra battery pouches and a power cable storage pouch



BACKPACK

ITEM

35-0037

SPECIFICATIONS

Dimensions: H 19 in. x W 9.5 in. x L 10.5 in.
Weight: 3.2 lbs.



APRU SOFT SHOULDER CARRIER



- Adds carrying capability by hand
- Provides additional impact protection
- Compatible with the Autonomous Portable Refrigeration Unit (APRU) - *sold separately*

SPECIFICATIONS

Dimensions: H 19.5 in. x W 10 in. x L 11 in.
Weight: 3 lbs.

SHOULDER CARRIER

ITEM

35-0038

DELTA ICE - 2L SMART BLOOD COOLER - 2 LITER



Contents not included

- Durable, rugged design, wipeable and easy to clean
- Insulation protected from drops and punctures
- Records & Stores Internal Temperature
- Displays Real-Time Temperature Reading
- Bluetooth and Cellular Connectivity. Alerts warn you of temperature excursions
- Aluminum insert featuring integrated, built-in NIST traceable sensors
- Maintains two 500 ml of whole blood for 48 hours at 35.6 - 42.8°F, and eliminates the need for glycol



SPECIFICATIONS

Dimensions: H 9.5 in. x W 9.7 in. x L 11.5 in.
Weight (unloaded): 10.5 lbs.

Insert Only:

Dimensions: H 5.75 in. x W 6.25 in. x D 8.24 in.
Weight: 5 lbs.

DELTA ICE

OPTION	ITEM #
Black	35-0040
Blue	35-0041
Red	35-0042
Insert Only	35-0043

Circulation & Vascular Access

LIFEWARMER™ Quantum™ IV Fluid & Blood Warming System

- **Lightweight** The Quantum™ has the lightest operational weight of any portable blood & fluid warming system (Total Operational Weight: 22 oz)
- **Powerful** Warms two units of 39°F (4°C) refrigerated whole blood at 100mL/min+ on a single charge
- **Compact** Fewest components to carry (three)
- **Fast** Warm-up time of only 24 seconds
- **Intuitive** The thermal transfusion/infusion set looks, handles, and sets up/primes like any conventional IV set
- **Safe** Tested safe for blood, blood products & IV solutions
- **Effective** The system is plug-and-play for ease of setup under stressful conditions such as tactical field care (only two connection points, with no plug orientation or set-up order requirements – system status constantly visible at a glance)
- **Complete** Does not require additional standard transfusion or infusion set – all other warmers do
- **Tactical** Only system FDA cleared allowing easy disable of audible alerts to protect operators in combat

SPECIFICATIONS

Total Operational Weight: 22 oz
(Quantum™ Controller, Battery & TTS-B)

Battery: 16 oz
Controller: 3 oz
Battery Charger: 22 oz
Controller Cover: 0.23 oz
Thermal Transfusion Set (TTS-B): 2.7 oz

COMPLETE QUANTUM IV & BLOOD WARMING SYSTEM CONTENTS:



QUANTUM COMPLETE SYSTEM

ITEM #	NSN #
35-0006	6515-01-693-7810

BLOOD or BLOOD PRODUCTS

From a **single battery charge**, capable of warming **1000 mL** of **39°F (4°C) Whole Blood** to an output temperature of **100.4°F (38°/- 2°C)**.
Warm up time: **24 seconds**.
Delivery Flow rate: **100mL/min**.

IV FLUIDS

From a **single battery charge**, capable of warming **1700 mL** of **68°F (20°C) IV fluids** to an output temperature of **100.4°F (38°/- 2°C)**.
Warm up time: **24 seconds**.
Delivery Flow rate: **200mL/min**.

TTS-B Thermal Transfusion Set

- Compatible with whole blood, component therapy (including FDP), colloids, and isotonics
- Similar size and weight to traditional IV tubing
- Intelligent sensing & heating elements integrated into the walls of the tubing set
- No disruption to the fluid path or laminar flow of the fluid
- Easy bi-directional connection to the Quantum Controller
- 20 gtt/mL drip rate

SPECIFICATIONS

TTS-B package (each): 6.5 in. x 11 in.
Tubing Length: 80 in.
TTS-B Wgt (each): 2.6 oz
Packaged Wgt (2 TTS-B & box): 11 oz



TTS-B (2 Pack)

ITEM #	NSN #
35-0005	6515-01-695-4444

TIS Thermal Infusion Set

- Similar size and weight to traditional IV tubing
- Intelligent sensing and heating elements integrated into the walls of the tubing set
- No disruption to the fluid path or laminar flow of the fluid
- Easy bi-directional connection to the Quantum Controller
- Compatible with standard IV fluids such as colloids and isotonics
- 20 gtt/mL drip rate

SPECIFICATIONS

TIS packaged (each): 6.5 in. x 11 in.
Tubing Length: 80 in.
Weight: TIS (each): 2.1 oz
Packaged (2 TIS & box): 10 oz



TIS (2 Pack)

ITEM #
35-0004

REV062121



Blood & Fluid Warming TRAINER SYSTEM

QUANTUM™ TRAINER SYSTEM & COMPONENTS

OPTION	ITEM #
Trainer System	35-0029
Trainer Tubing	35-0028
Trainer Controller	35-0027

- Allows medical trainers to instruct personnel on the proper use of the Quantum™ Blood and Fluid Warming System
- Designed for non-warming training purposes
- NOT FOR CLINICAL USE

CONTENTS:

- 2 x Trainer Tubing Sets
- 1 x Trainer Controller
- 1 x Battery
- 1 x BATTERY CHARGER

SPECIFICATIONS

Dimensions:

Trainer Tubing Package: W 6.5 in. x L 11 in. (ea.)
Tubing Length: 85 in.
Controller Body: W 1.25 in. x H 3 in. x D 0.75 in.
Controller Total Length: 51.25 in.

Weights:

Trainer Tubing Set: 2.7 oz.
Trainer Controller: 2.5 oz.
Trainer Controller (Pkgd): 5 oz.
Battery: 16 oz. | Charger: 22 oz.
Total Trainer Operational Weight: 22 oz.



> RESUSCITATE

500 ml in < 2 minutes
Stay with Your Patient

Deliver first bolus quickly to reverse hypotension in patients with shock

- Provide active resuscitation
- Maintain control over delivery

> RESTORE

**Improve Perfusion /
Reverse Shock**

Assess using vitals & bedside markers of perfusion, decide if additional fluid bolus is needed.

Early fluids more rapidly reverse shock & **improve outcomes**.¹⁻³

> REASSESS

Rapid Decision / Disposition

Each bolus delivered provides an opportunity to evaluate and determine next steps.

Faster resuscitation may lead to:

- Less mechanical ventilation⁴
- More rapid disposition
- **Less donning/doffing of PPE**

LifeFlow® PLUS BLOOD & FLUID INFUSER

- Hand-operated, rapid infusion device that allows health care providers to deliver fluid & blood products quickly and effectively
- Deliver 500ml of fluid in less than 2 minutes through a 20G IV catheter
- Deliver 1 unit of whole blood in 2 minutes
- Goes from the field to the hospital
- Can be operated with one hand
- Handle and tubing packaged together for quick and easy set-up
- Designed to protect the syringe from bacterial contamination
- Average time to set up and prime: less than 3 minutes

SPECIFICATIONS

12 in. x 7.5 in. x 2.75 in.
Weight: 0.95 lb

LIFEFLOW PLUS BLOOD & FLUID INFUSER

ITEM

Rx ONLY 30-0287


Replacement Tubing
(30-0288) available at
NARescue.com

¹ Lee, S. J., et al. (2014). "Increased fluid administration in the first three hours of sepsis resuscitation is associated with reduced mortality: a retrospective cohort study." Chest 146(4): 908-915.

² Williams, J. M., J. H. Greenstade, C. A. Dymond, K. Chu, A. F. T. Brown and J. Lipman (2018). "Characteristics, treatment and outcomes for all emergency department patients fulfilling criteria for septic shock: a prospective observational study." Eur J Emerg Med 25(2): 97-104.

³ Williams, J. M., J. H. Greenstade, C. A. Dymond, K. Chu, A. F. T. Brown and J. Lipman (2018). "Characteristics, treatment and outcomes for all emergency department patients fulfilling criteria for septic shock: a prospective observational study." Eur J Emerg Med 25(2): 97-104.

⁴ Leisman D, et al. (2017). Patterns and Outcomes Associated With Timeliness of Initial Crystalloid Resuscitations in a Prospective Sepsis and Septic Shock Cohort. Crit Care Med. Oct;45(10):1596-1606.

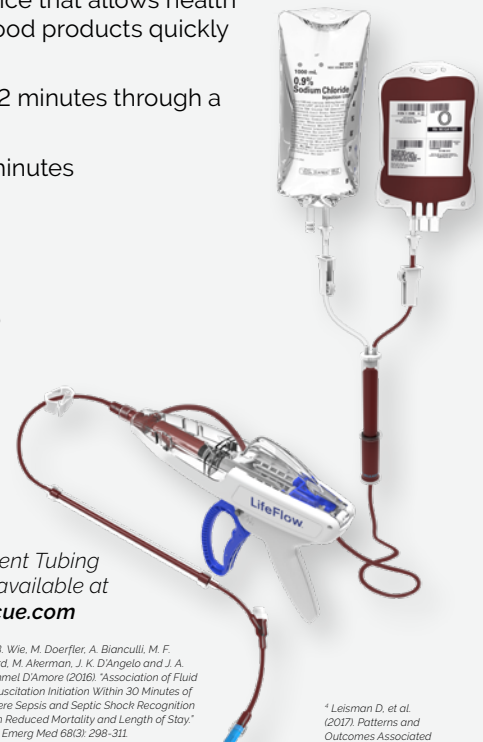



Hypothermia Prevention & Management

HPMK® Hypothermia Prevention & Management Kit

- Includes reinforced Heat Reflective Shell (HRS™) of 4-ply, composite fabric with a protected non-conductive reflective layer
- Built-in hood, integrated fluid absorption pad & top-to-bottom taper providing excellent thermal insulation
- Continuous 1.5 in. hook & loop closures provide rapid 360° access to the casualty with minimal exposure to the elements
- Self-heating, 4-cell liner (Ready-Heat™ Blanket) provides 10 hours of continuous dry heat up to 125°F (52°C)

SPECIFICATIONS

Shell Packaged:

H 6.75 in. x W 10.5 in. x D 5.5 in.

Shell Open: L 78 in. x W 43 in.

Weight: 3 lb 8 oz

US Patent: 7,766,950 B2

HPMK® CONTENTS:



HPMK

ITEM #	NSN #
80-0027	6515-01-532-8056

NOTE: Not for sale in Canada



NOTE: Litter is sold separately

HRS® Heat Reflective Shell

As an individual unit, reinforced 4-ply, composite shell with built-in hood, weight tested to support 250+ lb

SPECIFICATIONS

Shell Packaged: L 10 in. x W 6.75 in. x D 2.5 in.

Open: L 78 in. x W 43 in.

Weight: 1 lb 13 oz

US Patent: 7,766,950 B2

HRS - HEAT REFLECTIVE SHELL

ITEM #	NSN #
ZZ-0051	6532-01-627-7888



NAR® READY HEAT™ 4-CELL

- Blanket begins heating immediately upon opening, no activation required
- Approved for use on fixed wing and rotary aircraft by the US military
- Provides 10 hours of continuous dry heat up to 125°F (52°C)

SPECIFICATIONS

Packaged: L: 15 in. x W: 12.75 in. x D: 1 in.

Blanket: L: 26 in. x W: 22 in. x D: 0.5 in.

Weight: 188 lbs

NAR READY HEAT 4-CELL

ITEM

ZZ-1467



Hypothermia Prevention & Management



HPMK-I™ Hypothermia Prevention & Management Kit - Insulated

- Includes Heat Reflective Shell - Insulated (HRS-I™) constructed of 5-ply composite fabric with protected non-conductive reflective layer and NEW Climashield® insulation layer top & bottom
- Includes NAR Ready Heat® 4-Panel Blanket begins heating immediately upon opening

HPMK-I™ CONTENTS:



SPECIFICATIONS

Packaged - H 8 in. x W 11 in. x D 6 in.

Weight: 4 lbs 14 oz

Heat Reflective Shell-Insulated™ (HRS-I™)

(Open): L 78 in. x W 43 in.

U.S. Patent 7,766,950 B2, U.S. Patent: Pending

HPMK-I™

ITEM #	PRODUCT	NSN #
80-1192	HPMK-I™	---




NAR READY HEAT® 4-PANEL BLANKET



NEW NON-WOVEN MATERIAL BACKING



5-PLY COMPOSITE FABRIC WITH CLIMASHIELD® INSULATION

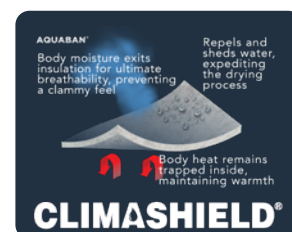

360° RAPID ACCESS HOOK & LOOP CLOSURES

HRS-I™ Heat Reflective Shell - Insulated™

- Insulation value exceeds that of the Army Green Sleeping Bag in service today (HRS-I offers total clo of 2.79 compared to Army Green Sleeping Bag total clo of 2.52)
- Highly compressible insulation has excellent rebound after opening and features Aquaban® technology making it hydrophobic
- Rapid 360° access through 1.5 in. continuous hook & loop closures located along the perimeter of the covering

HRS-I™

ITEM #	PRODUCT	NSN #
ZZ-1452	HRS-I™	---



SPECIFICATIONS

Packaged: H 11.25 in. x W 9 in. x D 4.25 in.

Size: L 78 in. x W 43 in.

Weight: 2.81 lbs

Customized Trauma Kits

Customizable Solutions built to meet YOUR requirements

If you haven't found exactly what you need, our team of **subject matter experts** at North American Rescue are standing by to guide you to the **ideal lifesaving solution** tailored to meet all your specifications.

Whether it's the essential gear needed for individuals or customized kits for multiple casualties, a **North American Rescue Specialty Kit** is your answer to the precise aid kit you and your organization need to **increase survivability** should injury occur.

High
Quality
Tactical
Nylon

-OR-

Vacuum
Bag

Call or Email us TODAY to
discover all your options for the
tailored lifesaving solution
YOU can live with!



Explore all the possibilities at
NORTH AMERICAN RESCUE®

www.NARescue.com • 888.689.6277

First responders have increasingly become **the target of violence** in recent history.

Nationwide, **violent crime is rising** at an alarming rate, translating to a higher call volume. With the **sharp decline** in the **police presence** which has historically provided scene safety for Fire/EMS, First Responders are experiencing a **dramatic increase** in **dangerous** and **potentially violent situations** in their everyday operations. Unfortunately rescuers in the Fire/EMS arena can be perceived as just another target for potential life threatening violence along with their law enforcement counterparts.

Now more than ever, your **personal safety is paramount** in the performance of your duties. Wearing **ballistic protective equipment (BPE)** is one of the best ways **to help protect yourself** in a violent environment.

PROTECT & SAVE Yourself Lives...

Law enforcement transitioned to wearing body armor on a daily basis in the 1980s because they realized it was an essential component of their daily equipment.

It saves lives! In the same way, BPE is critical for Rescue Task Force (RTF) and Tactical Emergency Medical Service (TEMS). Fire/EMS should now consider Ballistic Protection Equipment for daily wear on any call that could result in violence or an unsafe environment.

Protecting Yourself means you can Save others AND Live to serve another day!



◀ NIJ Level III + Enhanced Ballistic Rifle Plates*

NIJ ▶ 0101.06 Level IIIA Responder PPE Vest* Ballistic Protection



▲ Abdominal Extender & Custom ID Patches (each sold separately)



▲ Chest Pouch Medical Kit



Rescue Task Force Rapid Response Kit



▲ MCI Rapid Response Kit



▲ MCI WALK

*Ballistic protective equipment (BPE) products are controlled by the Export Administration Regulations (EAR). It may not be exported, sold, transferred, transported, or shipped outside of the U.S., either in its original form or after being incorporated into other items, without first obtaining approval from the U.S. government or as authorized by U.S. law and regulations.



Explore all your options at
NORTH AMERICAN RESCUE®
www.NARescue.com • 888.689.6277

www.NARescue.com • 888-689-6277 | **21**

RAPID RESPONSE KIT RTF Edition

- Equipped to provide initial triage of a mass casualty incident and provide treatment for 4 or more patients
- Includes 4 C-A-T® tourniquets with holders and nylon dispenser with 4 rolls of color-coded triage tape

RAPID RESPONSE -RTF

COLOR	ITEM #
Black	85-0716
Red	85-0717

SPECIFICATIONS

H 12 in. x W 12 in. x D 6 in.
Shoulder strap: 28 in. to 50 in.
Weight: 6 lb 4.5 oz

RAPID RESPONSE KIT - RTF EDITION:



Rapid Response Kits may be modified upon request. Contact your sales rep for info.

Triage Tape, Set of 4
(Item# ZZ-0591) also available

RAPID RESPONSE BAG RTF Edition

- Modular, personalized organization of mission essential gear (sold separately)
- Rugged nylon shoulder bag features adjustable, padded strap with quick release buckle
- MOLLE'd on the exterior are 4 nylon TQ Holders and 1 Triage Tape Dispenser (TQ's and Triage Tape sold separately)
- Multiple storage capabilities including loop-tape/MOLLE attachment system, 2 extra pockets & elastic storage loops with color-coded labels to accommodate up to 8 chem lights (not included)

SPECIFICATIONS

Bag Closed: H 12 in. x W 12 in. x D 6 in.
Adjustable shoulder strap: 28 in. to 50 in.
Weight: 2.05 lb

RAPID RESPONSE BAG - RTF

COLOR	ITEM #
Black	80-0412
Red	80-0413



Multiple Injury Trauma Kits

MCI-TKD

Mass Casualty Incident
Throw Kit Dispenser

- Contains an inner main compartment specifically designed for quick deployment of up to eight (8) point of injury kits
- One-handed quick release, silent zippers
- G-Hook backpack straps can switch from 2 shoulder straps to bandolier style.
- Easy backpack carry for mass casualty events
- Quickly deploy bleeding control kits in emergency situations
- Side pockets have C-A-T® Tourniquets and Emergency Trauma Dressings (ETD®) for maximum preparedness for multiple injuries
- Front pockets for Litter, Breathing, and wound packing scenarios



Every MCI-TKD includes these core items:



PLUS+

BASIC
also includes:

w/BCD
also includes:

w/CELOX™
also includes:



MCI-TKD Bleeding Control Throw Kit includes:



PLUS+

BASIC
also includes:

w/BCD
also includes:

w/CELOX™
also includes:



SPECIFICATIONS

H 21 in. x W 13 in. x D 8 in.
Weight: 15 lbs.

MCI-TKD

COLOR	LEVEL	ITEM #
● BLK	Basic	80-1749
● RED	Basic	80-1750
● BLK	w/ BCD	80-1751
● RED	w/ BCD	80-1752
● BLK	w/ CELOX™	80-1753
● RED	w/ CELOX™	80-1754



NEW ITEM!

Multiple Injury Trauma Kits

MINI TRAUMA SHEET BAG

Designed for portability, visibility, and easy access

Introducing the NAR® Mini Trauma Sheet, a portable, compact carry sheet that can carry life saving equipment. Contents can be customized to meet your mission requirements. The Mini Trauma Sheet can be rolled up to be carried to an incident, or hung in an ambulance, vehicle, or aircraft for easy access to medical supplies.

The Mini Trauma Sheet features clear vinyl pouches for easy identification of contents, and each pouch can be accessed through the large zipper opening. Pouches have loop fields enabling the use of hook and loop name tape labeling such as Bleed, Breathing, MISC, etc (name tapes are not included).



SPECIFICATIONS

Closed: 13 in. W x 9.75 in. L x 2.75 in. D

Open: 13 in. W x 19 in. L x 1 in. D

Weight: .78 lbs.



Full Size

A compact version of NAR's Medic/Trauma Sheet



MINI TRAUMA SHEET KIT



- Can be rolled up and carried with the attached handles, or opened and hung using the sewn in D-Rings
- Four large pouches and one lower profile zippered pocket for maximum capacity
- Compact design for smaller vehicles, ambulances, etc.
- Ability to pack out however you like for best possible outcomes
- Hook and Loop attachment for name tape additions for better quality quick visual access

MINI TRAMA SHEET

COLOR	Option	ITEM #
BLK	Basic	80-1808
COY	Basic	80-1809
RED	Basic	80-1807
BLK	w/ BCD	80-1805
COY	w/ BCD	80-1806
RED	w/ BCD	80-1804
BLK	Bag Only	80-1721T
COY	Bag Only	80-1722T
RED	Bag Only	80-1723T

SPECIFICATIONS

Closed: W 13 in. x L 9.75 in. x D 5 in.
Open: W 13 in. x L 19 in. x D 2 in.
Weight: 2.13 lbs.

PLUS+
with
Bleeding Control Dressing
also includes:



* RED Mini Trauma Sheet Kits are packed with
Orange C-A-T® Tourniquets and RED/BLK name patches.

Multiple Injury Trauma Kits

Mass Casualty Incident **MCI WALK**

- All-in-one trauma treatment & evacuation platform for treatment of multiple trauma patients
- Choose from two different litter configurations: 8 non-rigid, disposable QuikLitters or 1 full-size, quad-fold, lightweight MedEVAC Litter
- Pocket system holds self- or buddy-aid point-of-wounding medical supplies
- Bag features back pack straps for easy transport from the vehicle to the point of need
- Ideal for use as a vehicle trauma kit

(80-0565/80-0458) Red MCI WALKs contain orange C-A-T® Tourniquets

(#80-0458) Red MCI WALK with MedEVAC Litter contains an orange MedEVAC Litter

MCI WALK

BAG COLOR	LITTER OPTION	C-A-T	ITEM #
● Black	8 QuikLitters	Black	80-0564
● Red	8 QuikLitters	Orange	80-0565
● Black	MedEVAC Litter (Black)	Black	80-0457
● Red	MedEVAC Litter (Orange)	Orange	80-0458

SPECIFICATIONS

H 22.75 in. x W 16.5 in. x D 11.5 in.

Weight with 8 QuikLitters: 18 lb 12 oz

Weight with MedEVAC Litter: 29 lb 15 oz

MCI WALK with QUIKLITTERS:

 QTY 8 C-A-T® Tourniquet	 QTY 8 ETD® 6 in. Dressing	 QTY 5 Combat Gauze® LE	 QTY 8 NAR® Gauze, Compressed	 QTY 5 Pre-lubed NAR® Naso, 28F
 QTY 8 HyFin® Vent Chest Seal, Twin Pk	 QTY 8 NAR® ELITE Litter	 PAIR 8 Nitrile Gloves, Lg	 QTY 20 TCCC Triage Card	 QTY 5 Permanent Marker, Large

MCI WALK with MEDEVAC LITTER:

 QTY 8 C-A-T® Tourniquet	 QTY 8 ETD® 6 in. Dressing	 QTY 5 Combat Gauze® LE	 QTY 8 NAR® Gauze, Compressed	 QTY 5 Pre-lubed NAR® Naso, 28F
 QTY 8 HyFin® Vent Chest Seal, Twin Pk	 QTY 1 MedEVAC Quad-Fold Litter	 PAIR 8 Nitrile Gloves, Lg	 QTY 20 TCCC Triage Card	 QTY 5 Permanent Marker, Large

Bag only option also available at **NARescue.com**



online video available

< scan to view

SLISHMAN TRACTION SPLINT



STS-COMPACT



STS-GEN 2

- One of the smallest and lightest traction splints on the market
- One size fits all, from small children to large adults
- Can be applied in under a minute with minimal training

- Single connected device—no loose parts, which may get lost
- Application possible even in cases of lower leg/foot amputation
- Radiolucent

SLISHMAN TRACTION SPLINT

OPTION	ITEM #	NSN #
STS-Gen 2	50-0047	----
STS-Compact	50-1008	6515-01-694-7603

SPECIFICATIONS

Gen 2: Packaged: L 24 in. x W 4.75 in. x D 2.25 in.
Fully Extended: L 22.25 in. x W 4.25 in. x D 2.25 in.
Weight: 1.19 lbs.

Compact: Packaged: L 13 in. x W 3 in. x D 3 in.
Weight: 0.81 lbs.

Trusted by major EMS agencies including Boston EMS, Dallas Fire-Rescue, LA County Fire, and FDNY

PELVIC BINDER from CRO Medical

- BOA Precision Dial System
- Hook and loop retention for Pressure Delivery Device
- Approximately half the size in cubic inches of other binders
- Lightweight 4 in. double-back wrap
- Easy, one-hand tightening
- Easily achieves at least 33 lb (150N) of circumferential force
- Easily reduces and stabilizes open-book pelvic fractures using its patented tightening system



PELVIC BINDER

ITEM #
50-0044

SPECIFICATIONS

Standard Waist Size: 32 in. - 50 in. hips
Weight: 0.75 lb

SAM® PELVIC SLING for Pelvic Ring Fractures

- Force-controlled, ideal compression every time without the guesswork
- Rapid, simple, three-step application
- Radiolucent and MRI safe



SAM PELVIC SLING

ITEM #	NSN #
50-0004	6515-01-509-6866

SPECIFICATIONS

H 7.25 in. x W 7.75 in. x D 4 in.
Hip Circumference: 32 in. to 50 in.
Weight: 9.8 oz

NAR® ELITE™ - LITE Responder Litter

- Multi-purpose responder litter
- Perfect balance of size, weight, cost, and capability
- Features ten large carry handles with a load capacity tested to 750 lb
- Vacuum-packed to conserve space and fit in a cargo pocket
- Latex-free & water-resistant



SPECIFICATIONS

Payload 750 lb

NAR ELITE - LITE

ITEM #
60-0071

Litter Bed:
L 74 in. x W 28 in.
Weight: 13 oz

Heavy Duty Trauma **XSHEAR®**

- Hardened stainless steel blades that are twice as thick as most trauma shears, sharpened to a razor edge
- Patented blunt tip and curved design for gentler edge near the skin
- Soft touch inner handles for comfort and slip resistance
- Heavy duty center bolt that won't loosen over time

SPECIFICATIONS

L 7.5 in x 3.5 in. x 0.118 in.
Cutting area length: 2 in.
Weight: 4 oz

XSHEAR®

HANDLE	BLADE	ITEM #
Black	Black	70-0290
Blue	Silver	70-0291
Orange	Black	70-0292
Pink	Black	70-0293
Red	Black	70-0294
Coyote	Black	70-0415
Green	Silver	70-0424

Rigid **XSHEAR®** HOLSTER



Shears & Holster are sold separately

- Custom designed for the XSHEAR® Heavy Duty Trauma Shears
- MOLLE Compatible adjustable clip
- Flex lock design securely holds shears in place with a positive click
- Universal design allows shears to be inserted either direction
- Rugged injection molded ABS Plastic

RIGID XSHEAR® HOLSTER

ITEM

70-0295

SPECIFICATIONS

L 5 in. x W 3.75 in. x D 1.75 in. | Weight: 3.2 oz



Shears & Holsters sold separately

XSHEAR® Gen 2 **TACTICAL HOLSTER**



- Rugged molded ABS construction
- Flex-lock Dual Band retention system
- Smooth, rapid deployment
- Configurable for left or right carry
- Universal Belt Clip Bracket
- Discreet Carry Concepts clip included
- Made in the USA

SPECIFICATIONS

Box - L: 6 in. x W 4.25 in.
Weight: 5.76 oz.

XSHEAR® GEN 2 HOLSTER

ITEM

70-0425

Soft **XSHEAR®** HOLSTER

- Designed for XSshears but compatible with most shears
- Durable Spandex over soft foam with nylon interior
- Powder coated heavy steel belt clip
- Allows Left or Right hand carriage

SOFT XSHEAR® HOLSTER

ITEM

70-0314

SPECS

L 7 in. x W 3.5 in. x D 0.5 in.
Weight: 3.2 oz



SPECIFICATIONS

Shear Length: 19 in.
Closed Length: 5 in.
Weight: 5.8 oz

Leatherman® **RAPTOR® MULTI-TOOL**

- Integrated multi-tool specific to emergency medical situations
- High quality 420HC stainless steel medical shears
- Foldable shears come apart for easy cleaning
- Injection-molded polymer Holster on a rotating clip; MOLLE-compatible

1. 420HC Stainless Steel Folding Medical Shears
2. Strap Cutter
3. Oxygen Tank Wrench
4. Ring Cutter
5. Carbide Glass Breaker
6. 5 cm Ruler
7. Holster with Clip

RAPTOR® MULTI-TOOL

HANDLE	ITEM #	NSN #
Black	70-0014	5110-01-626-8766
Coyote	70-0265	----
Org/Blk	70-0327	----
Blu/Blk	70-0264	----
Red/Blk	70-0266	----
Blu/Ylw	70-0326	----

ARMADILLO

Medication & Storage Case

- Compact solution for keeping critical medications safe, secure and easily accessible
- Features Automatic Pressure Equalization Valve
- Hardware constructed of Marine grade 316 SS
- Case is water-resistant and crush-proof
- Hook & loop attachment system allows storage compartment to be tailored to hundreds of custom configurations

SPECIFICATIONS

Ext: L 7.95 in. x W 5.1 in. x D 2.14 in.
Int: L 6.45 in. x W 3.87 in. x D 1.75 in.
Weight: 10.5 oz

ARMADILLO MEDICATION CASE * NATO NSN

LID STYLE	ITEM #	NSN #
● Black	ZZ-0081	8145-01-573-2533
○ Clear	ZZ-0259	6545-01-594-9399*



Clear Lid
ZZ-0259 only
includes Loop Liner
on the bottom

also available: Half Inch Hook Discs with pressure sensitive adhesive
Set of 250 • Item# ZZ-0258 • NSN# 8145-01-617-9075*

SCISSOR LEASH™

- Salt water resistant
- Self cleaning
- Helps keep track of your shears

SPECS

H 6 in. x W 16 in. x D 0.6 in.
Cable Length - 36 in.
Weight: 16 oz

SCISSOR LEASH

ITEM #	NSN #
70-0005	6515-01-540-7226



6 oz MOLLE Retractable GEAR LEASH

- MOLLE Mount compatible with 360° rotation
- High break strength / Salt water resistant
- 36 in. Extension / Spectra Nylon Line
- Self-Flushing
- Stainless Steel Spring & Hardware

GEAR LEASH

ITEM #
70-0155



SPECS

Retracted:
H 3.5 in. x W 1.6 in. x D 1.0 in.
Cable Length - 36 in.
Weight: 2.0 oz

Blue RESPONDER NITRILE GLOVES

- Form-fitting, ambidextrous design with textured fingertips
- Thick palm ensures durability
- Textured for improved wet and dry gripping
- Longer and thicker than standard nitrile gloves
- Available in pre-packaged pairs (kits) or in bulk (box)

SPECIFICATIONS

Bag (25 pair): L 4.75 in. x W 6 in. x D 4 in.
Weight: 1 lb 1 oz

Box (Bulk): L 11 in. x W 5 in. x D 3.75 in.
Weight: 2 lb 4 oz

Blue Responder BULK

SIZE	PKG QTY	ITEM #	SIZE	PKG QTY	ITEM #
S	Box 50 pair	70-0315	L	Box 50 pair	70-0317
M	Box 50 pair	70-0316	XL	Box 50 pair	70-0318



COLOR:
Blue



Blue Responder PAIRS

SIZE	PKG QTY	ITEM #
S	Bag 25 pair	ZZ-1059
M	Bag 25 pair	ZZ-1060
L	Bag 25 pair	ZZ-1043
XL	Bag 25 pair	ZZ-1061

* also available in black at NARescue.com

Take your training to the
NEXT LEVEL
with superior **REALISM**

Tactical
Operation
Medical Manikin
SERIES

TOMManikin®

BASIC | BLAST | GSW | WATER | CONFINED SPACE | BURN | CBRNE | FEMALE

Interactive Animatronics:

- Lifelike Breathing Mechanisms
- Functional Anatomical Airway
- Variable Radial Pulse Rates & Pressure
- High Fidelity Bleeding from Multiple Wounds
- Customizable, Multi-language Voice Tracks
- Rugged Softech® Skin in a variety of skin tones

Wireless Controls:

- MILSPEC Wireless Remote Instructor Tablet & Software
- Stand alone point-to-point communications with .5 mile line of site range

Flexible Wound Configurations:

Basic | Blast | GSW | Water | Confined Space | Burn | CBRNE
Custom configurations including Female Anatomy on request

Medical Procedural Capabilities:

MASSIVE HEMORRHAGE:

- Tourniquet Application (proximal, junctional & abdominal)
- Active Bleeding Control by Direct Pressure on Pressure Points
- Wound Packing

ADVANCED AIRWAY:

- Nasopharyngeal/Oral/Supraglottic Airway Insertion
- Endotracheal Intubation
- Surgical Airway and Tissue Dissection with Realistic Bleeding

RESPIRATIONS:

- Chest Seal Application
- Bilateral Needle Decompression
- Chest Tube Capable

CIRCULATION:

- IO (sternal & humoral)
- Antecubital IV Access on uninjured upper limbs



NORTH AMERICAN RESCUE®

www.NARescue.com • 888.689.6277

For more info, contact **NAR Simulation Sales** | manikins@narescue.com | **866.839.3455**

CIVILIAN EMS COURSE OFFERINGS

ENROLL TODAY!



IMMERSIVE



HANDS-ON



INSTRUCTOR
LED

NAR C3

42-hour in-house course for IBSC critical care course certification or re-certification (FP-C or CCP-C) and national registry re-certification. Students will learn or refresh all skills and knowledge needed to certify/re-certify, and will earn 42 hours of IBSC and National Registry Continuing Education.

PHTLS

The Pre-Hospital Trauma Life Support (PHTLS) Provider Course was developed by NAEMT for EMTs, AEMTs, Paramedics, Nurses, Physicians Assistants, and Physician. This is a 2-day in-house course based on a philosophy stressing the treatment of the multi-system trauma patient as a unique entity with specific needs. PHTLS promotes critical thinking as the foundation for providing quality care through hands on skills and high-quality lectures.

Enhanced HEMCON

8-hour course oriented towards non-medical members of the community. Students are trained to identify and treat massive hemorrhage through lectures, hands-on skills, and realistic trauma scenarios. Students receive critical knowledge on how to respond to preventable bleeding deaths prior to arrival of EMS.

Initial EMT/AEMT Hybrid Course

NAR Training offers students entry level education and training for emergency medical professionals. Self-paced online learning combined with a hands-on skills week with clinical labs are designed to prepare students to pass the NREMT test while allowing for schedule flexibility. Completion of this course allows student to take National Registry exam.

EMT/AEMT/Paramedic Online Refresher Course

NAR Training offers students ongoing education and retraining for emergency medical professionals seeking to maintain and improve their practice competencies and stay ahead of new industry developments and technologies. Our online courses offer the flexibility vital for busy professionals in pursuit of career advancement. This training is CAPCE accredited and meets or exceeds the NREMT's requirements for Continuing Education.



NORTH AMERICAN RESCUE®

Education and Training

ACCREDITATIONS:



AUTHORIZED
TRAINING

AUTHORIZED
TRAINING
SITE



American
Heart
Association®

Now SAME DAY Cards Issued At End Of Course!

Discover more at
NAReducation.com



Course Information & Registration:

NAReducation.com
888-374-7808, Ext 101
registration@NARescue.com

MASSIVE HEMORRHAGE | AIRWAY | RESPIRATION | CIRCULATION | HYPOTHERMIA



NORTH AMERICAN RESCUE®

www.NARescue.com • 888.689.6277

Sign up to get updates and
special offers when you...

#JOIN THE MISSION

